



DICOM and Slicer: A Tutorial

Sonia Pujol, Ph.D.


Assistant Professor of Radiology

Director of 3D Slicer Training & Education

Brigham and Women's Hospital

Harvard Medical School

spujol@bwh.harvard.edu

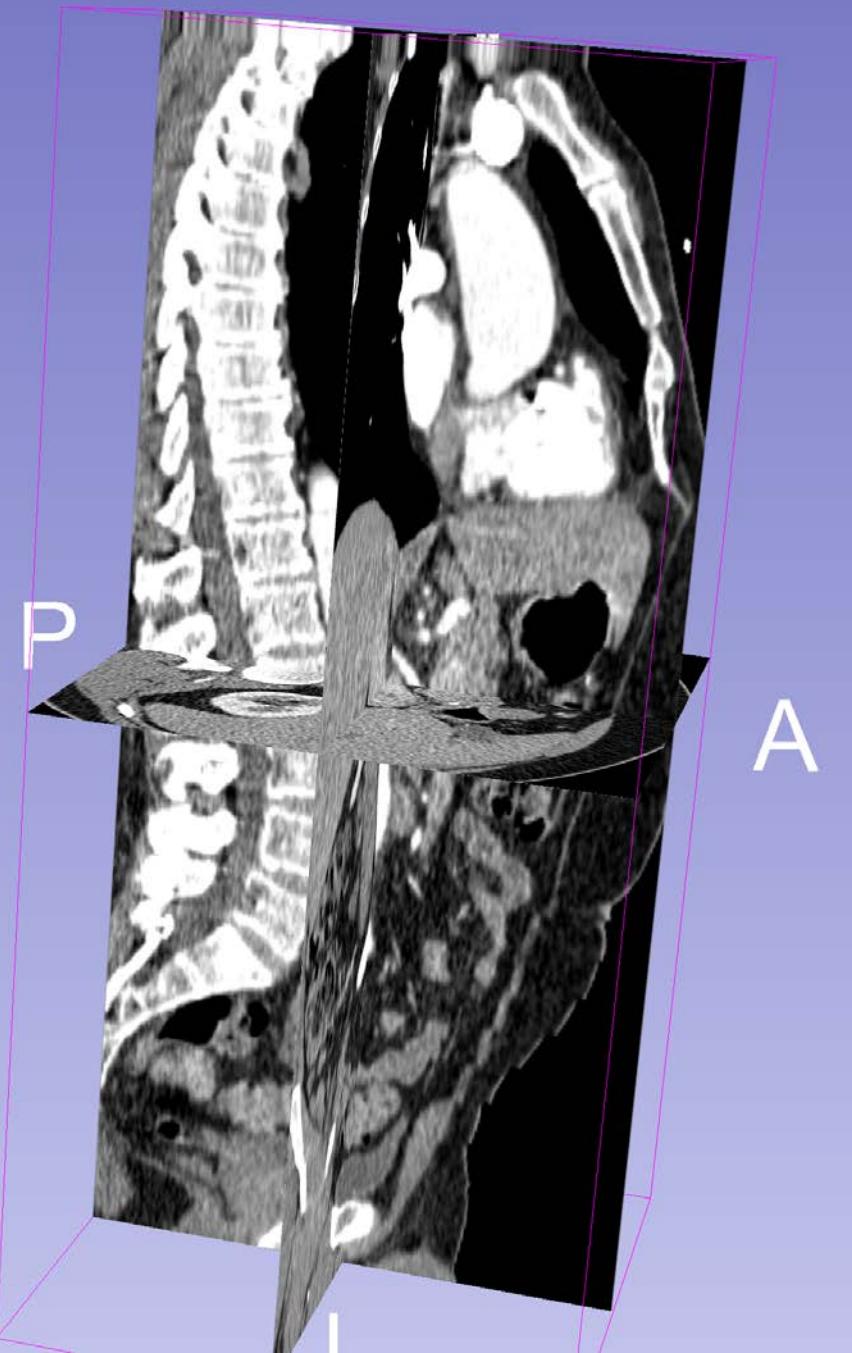

Objective

This tutorial provides a basic introduction to the DICOM standard, and shows how to visualize DICOM images in 3D Slicer version 5.0



Tutorial material



3D Slicer version 5.0

<https://download.slicer.org>

SlicerDICOMTutorialData

DICOM Torso CT

DICOM Breast MRI

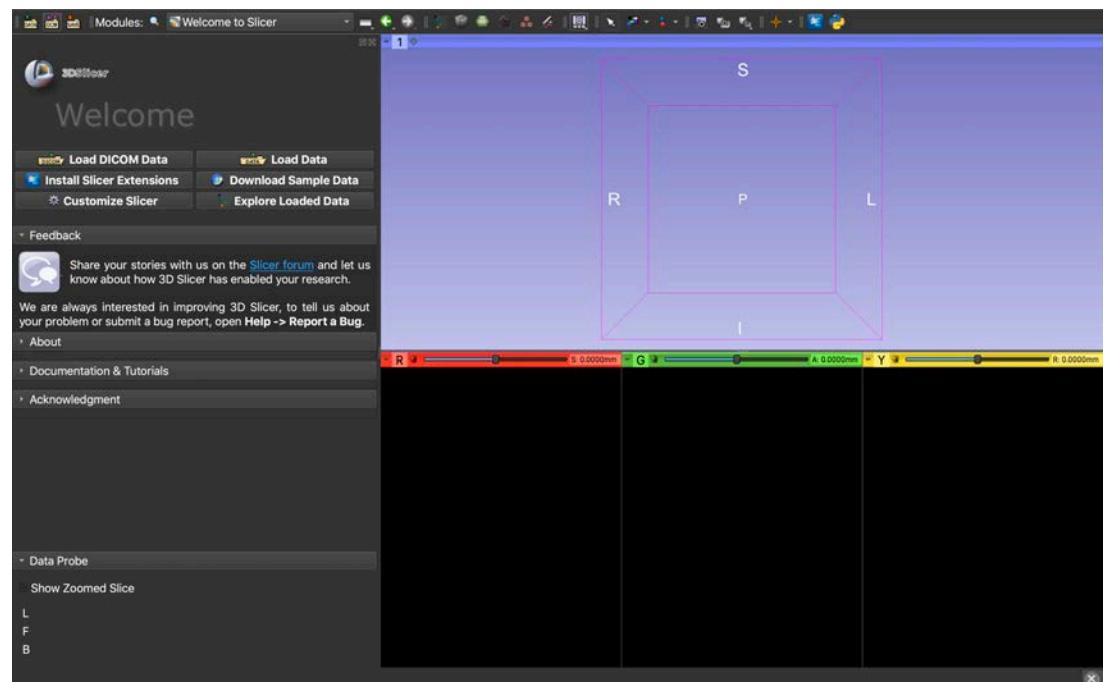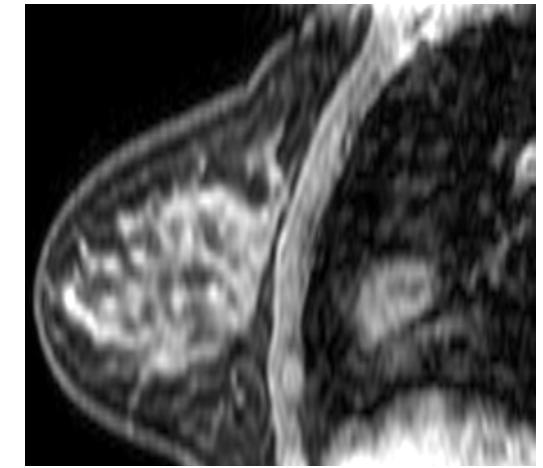

Tutorial Outline



Part 1: Introduction to DICOM



Part 2: Loading and Visualizing DICOM data in 3D Slicer version 5.0

Disclaimer

- 3D Slicer is a free open source software application distributed under a BSD style license
- The software is not FDA approved or CE-Marked, and is for research use only



Digital Imaging and Communications in Medicine


Part 1: Introduction to DICOM

Reproducible Science

- Reproducible science is critical to drive research and accelerate discoveries
- Open-source software tools such as 3D Slicer and data standards such as DICOM contribute to the reproducibility of scientific results in biomedical research



F.A.I.R. Principles

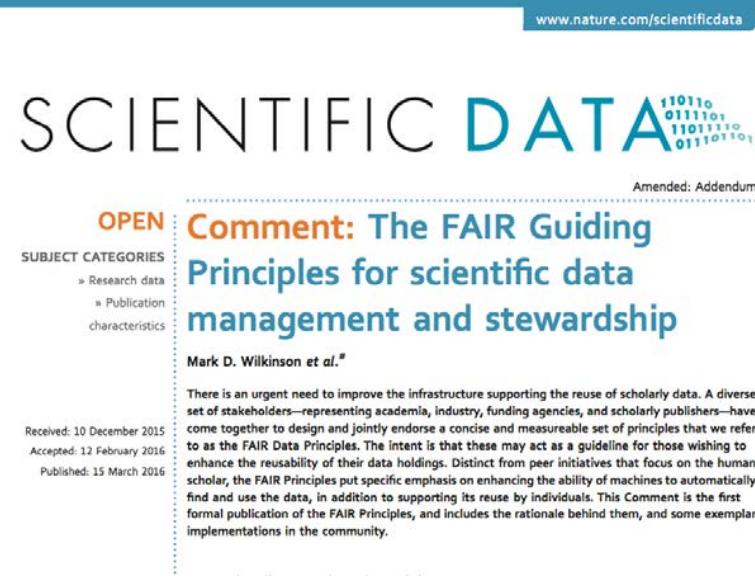

The FAIR Guiding Principles for scientific data management and stewardship.
Wilkinson et al. Sci. Data 2016
<http://go-fair.org/fair-principles>

- **Findable:** Data are easily findable
- **Accessible:** Users know how to access the data, including authentication and authorization
- **Interoperable:** Data can be integrated with other data and can interoperate with applications for storage and analysis
- **Reusable:** Data can be replicated or combined for new research

The DICOM standard

- DICOM (Digital Imaging and Communications in Medicine) is the international standard for handling, storing, printing and transmitting medical imaging data
- Clinical imaging equipment (CT scanners, MR scanners, X-Ray and ultrasound machines) generate DICOM files



DICOM History

- 1982: The American College of Radiology (ACR) and the National Electrical Manufacturers Association (NEMA) initiate standards for the interconnection of medical imaging devices
- 1985: Publication of the ACR-NEMA Digital Imaging and Communications Standards version 1.0
- 1988: Publication of the ACR-NEMA Digital Imaging and Communications Standards version 2.0
- 1993: Publication of the ACR-NEMA Standards version 3.0 also referred to as the Digital Imaging and Communication in Medicine (DICOM) standard

DICOM Today

- The DICOM standard is continuously being refined to address new community needs through multiple releases every year
- As of 2020/07/06, the DICOM standard is DICOM PS3 2020c and contains 4,000 pages
- DICOM Working Groups are established by the DICOM Standard Committee to develop the Standard for a particular modality, clinical domain or technical area (e.g. WG-16 Magnetic Resonance)

<http://dicomstandard.org>

FAIR Data and the DICOM Standard

www.nature.com/scientificdata

SCIENTIFIC DATA 

Amended: Addendum

OPEN

SUBJECT CATEGORIES
» Research data
» Publication characteristics

Comment: The FAIR Guiding Principles for scientific data management and stewardship

Mark D. Wilkinson et al.*


Received: 10 December 2015
Accepted: 12 February 2016
Published: 15 March 2016

There is an urgent need to improve the infrastructure supporting the reuse of scholarly data. A diverse set of stakeholders—representing academia, industry, funding agencies, and scholarly publishers—have come together to design and jointly endorse a concise and measurable set of principles that we refer to as the FAIR Data Principles. The intent is that these may act as a guideline for those wishing to enhance the reusability of their data holdings. Distinct from peer initiatives that focus on the human scholar, the FAIR Principles put specific emphasis on enhancing the ability of machines to automatically find and use the data, in addition to supporting its reuse by individuals. This Comment is the first formal publication of the FAIR Principles, and includes the rationale behind them, and some exemplar implementations in the community.



The DICOM Standard facilitates compliance of imaging data with FAIR principles

DICOM Data Model

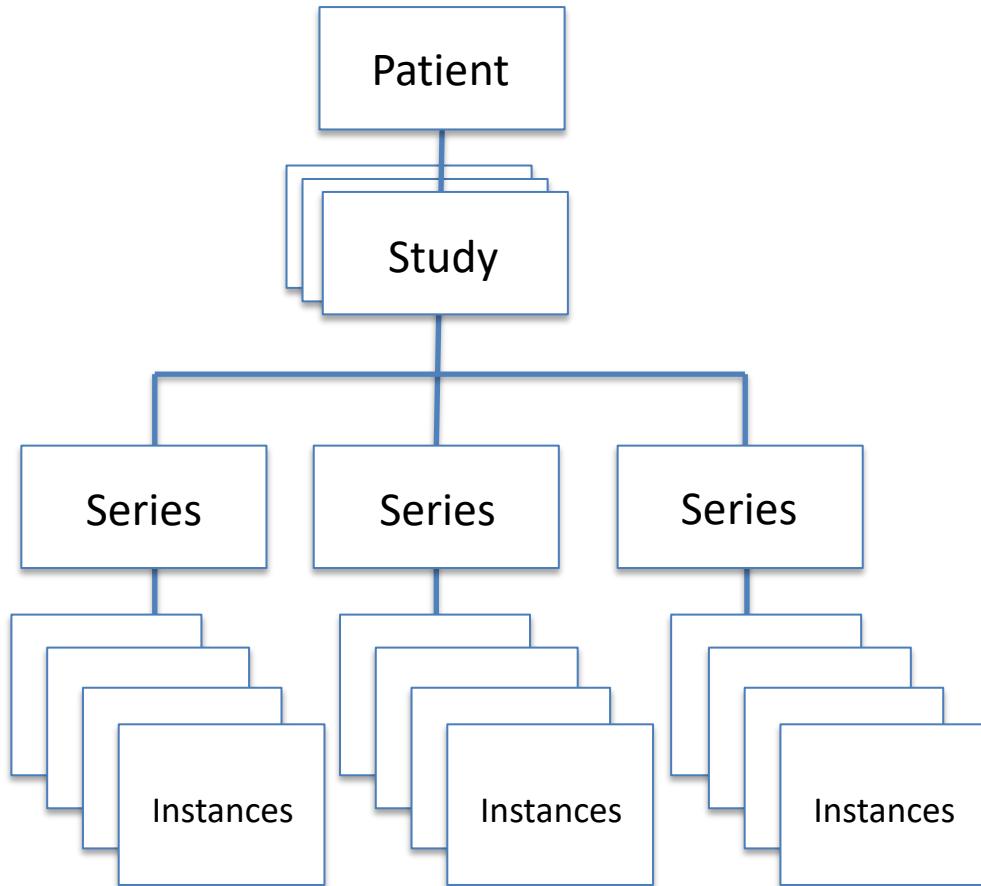

- In the DICOM data model, a **DICOM Study** consists of several **DICOM Series**, and each **DICOM Series** contains **DICOM Instances**
- Each of the **DICOM Studies**, **Series** and **Instances** are assigned a **Unique Identifier (UID)**

Example of DICOM instances: DICOM MRI Image data

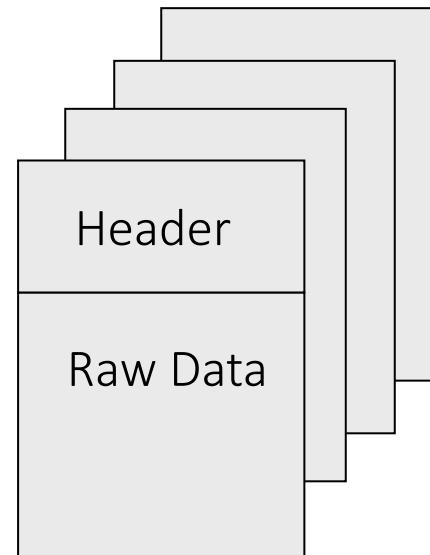
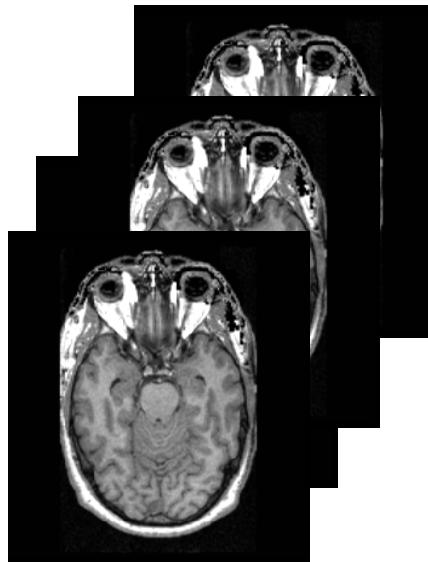


Image001.dcm

Image002.dcm

Image003.dcm


....

An MRI image is an example of DICOM Instance that consists of a **DICOM header** and an **image dataset**

```
0002,0000,File Meta Elements Group Len=148  
0002,0001,File Meta Info Version=256  
0002,0002,Media Storage SOP Class UID=1.2.840.10008.5.1.4.1.1.4.  
0002,0003,Media Storage SOP Inst UID=0.0.0.0.  
0002,0010,Transfer Syntax UID=1.2.840.10008.1.2.1.  
...  
0008,0060,Modality=MR  
0008,0070,Manufacturer=GE MEDICAL SYSTEMS  
0008,0080,Institution Name=1852796513  
0008,0081,City Name=1852796513  
0008,0090,Referring Physician's Name=1852796513  
0008,0092,?=1852796513  
0008,0201,?=-0500  
0008,1010,Station Name=1852796513  
0008,1030,Study Description=anon  
0008,103E,Series Description=anon  
0008,1040,Institutional Dept. Name=1852796513  
0008,1050,Performing Physician's Name=1852796513  
0008,1060,Name Phys(s) Read Study=1852796513  
0008,1070,Operator's Name=anon  
0008,1080,Admitting Diagnosis Description=1852796513  
0008,1090,Manufacturer's Model Name=GENESIS.SIGMA .....  
0010,0010,Patient's Name=anon  
0010,0020,Patient ID=anon  
0010,0030,Patient Date of Birth=00000000  
0010,0032,Patient Birth Time=000000  
0010,0040,Patient Sex=O  
0010,1010,Patient Age=000Y  
.....  
0028,0010,Rows=256  
0028,0011,Columns=256  
0028,0030,Pixel Spacing=0.937500 0.937500  
0028,0100,Bits Allocated=16  
0028,0101,Bits Stored=16  
0028,0102,High Bit=15  
0028,0103,Pixel Representation=1  
.....  
7FE0,0010,Pixel Data=131072
```

Example of DICOM header content

- The DICOM header contains metadata information on the patient, study and imaging data.
- Original data from the scanner tell users important elements about the acquisition
- The metadata information is accessible through DICOM tags
- DICOM tags uniquely identify DICOM attributes

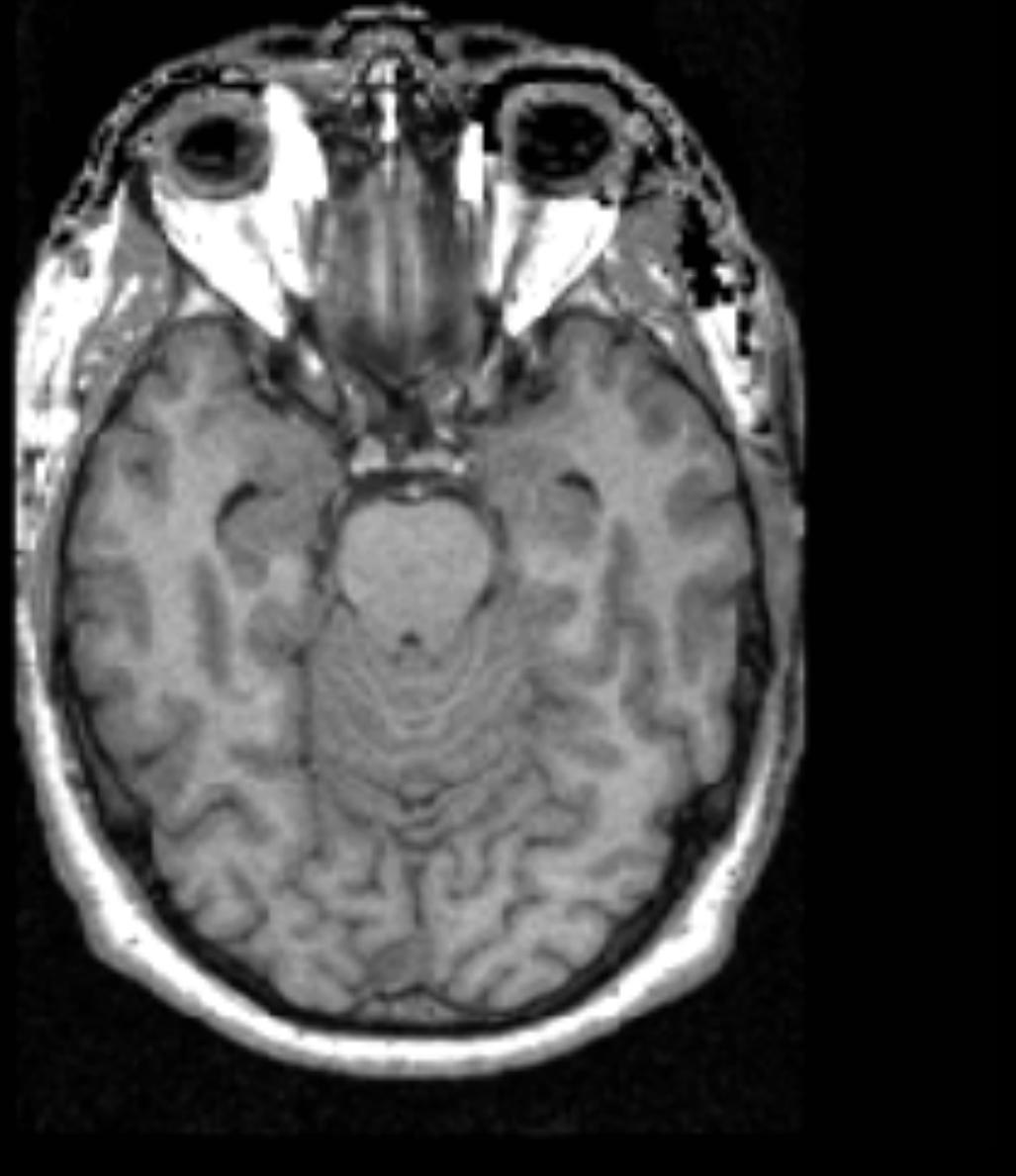

The Life Cycle of Medical Imaging Data - Sonia Pujol, Ph.D.
Functional Magnetic Resonance Imaging: Data Acquisition and Analysis
Harvard-MIT Health Sciences and Technology


0002
0002
0002
0002
0002

Physician and Study Information

0008,0060,Modality=MR
0008,0070,Manufacturer=GE MEDICAL SYSTEMS
0008,0080,Institution Name=1852796513
0008,0081,City Name=1852796513
0008,0090,Referring Physician's Name=1852796513
0008,0092,?=1852796513
0008,0201,?=0500
0008,1010,Station Name=1852796513
0008,1030,Study Description=anon
0008,103E,Series Description=anon
0008,1040,Institutional Dept. Name=1852796513
0008,1050,Performing Physician's Name=1852796513
0008,1060,Name Phys(s) Read Study=1852796513
0008,1070,Operator's Name=anon
0008,1080,Admitting Diagnosis Description=1852796513
0008,1090,Manufacturer's Model Name=GENESIS.SIGNA

0010,0010,Patient's Name=anon
0010,0020,Patient ID=anon
0010,0030,Patient Date of Birth=00000000
0010,0032,Patient Birth Time=000000
0010,0040,Patient Sex=O
0010,1010,Patient Age=000Y
.....

0028,0010,Rows=256
0028,0011,Columns=256
0028,0030,Pixel Spacing=0.937500 0.937500
0028,0100,Bits Allocated=16
0028,0101,Bits Stored=16
0028,0102,High Bit=15
0028,0103,Pixel Representation=1
.....
7FE0,0010,Pixel Data=131072

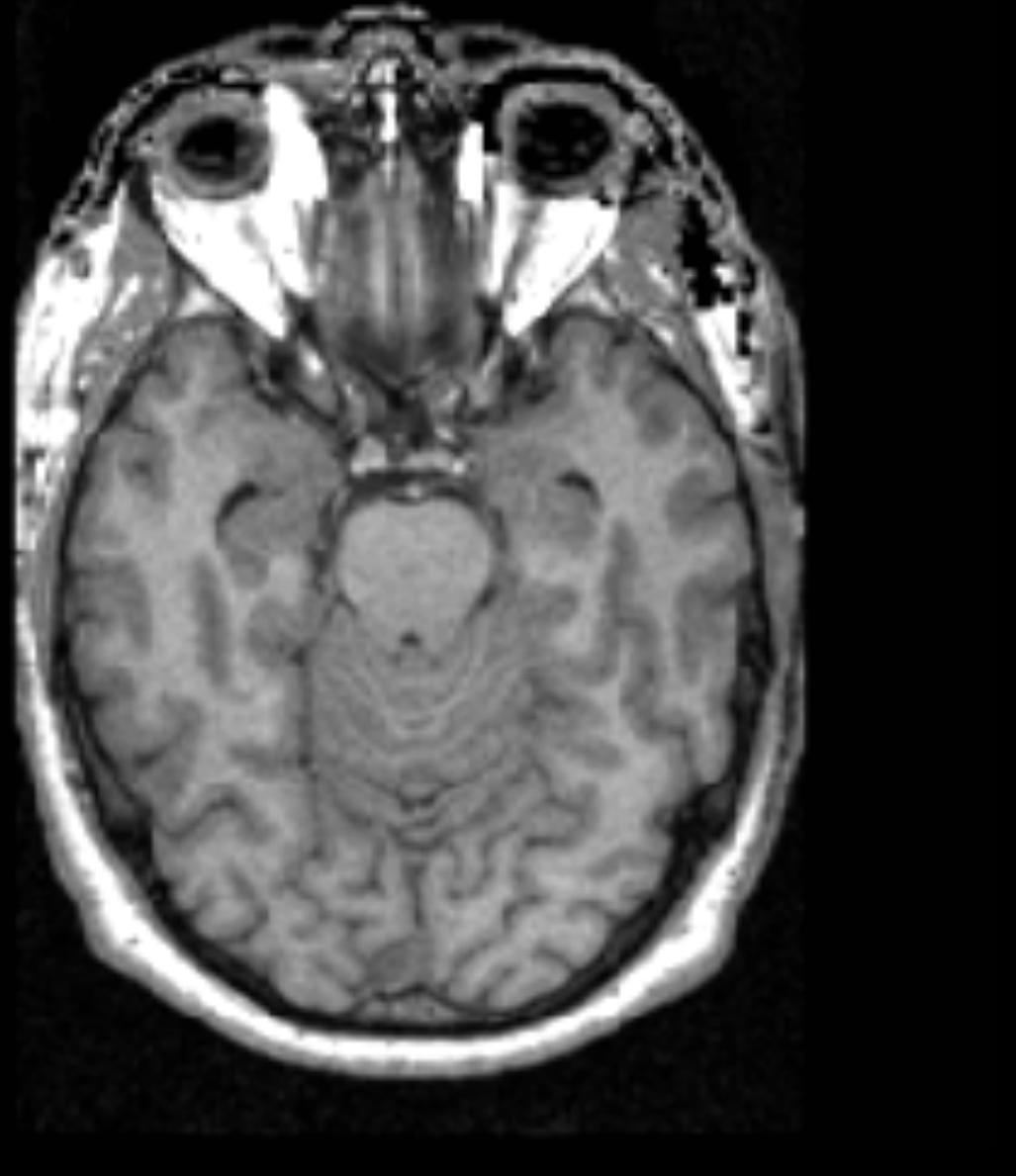

0002,0000,File Meta Elements Group Len=148
0002,0001,File Meta Info Version=256
0002,0002,Media Storage SOP Class UID=1.2.840.10008.5.1.4.1.1.4.
0002,0003,Media Storage SOP Inst UID=0.0.0.0.
0002,0010,Transfer Syntax UID=1.2.840.10008.1.2.1.

...
0008,0060,Modality=MR
0008,0070,Manufacturer=GE MEDICAL SYSTEMS
0008,0080,Institution Name=1852796513
0008,0081,City Name=1852796513
0008,0090,Referring Physician's Name=1852796513
0008,0092,?=1852796513
0008,0201,?=-0500
0008,1010,Station Name=1852796513
0008,1030,Study Description=anon
0008,1035,Series Description=
0008,
0008,
0008,
0008,
0008,1000,Admitting Diagnosis Description=1852796513
0009,1000,Manufacturer's Model Name=GENEPI SIGNA

0010,0010,Patient's Name=anon
0010,0020,Patient ID=anon
0010,0030,Patient Date of Birth=00000000
0010,0032,Patient Birth Time=000000
0010,0040,Patient Sex=O
0010,1010,Patient Age=000Y

.....
0028,0010,Rows=256
0028,0011,Columns=256
0028,0030,Pixel Spacing=0.937500 0.937500
0028,0100,Bits Allocated=16
0028,0101,Bits Stored=16
0028,0102,High Bit=15
0028,0103,Pixel Representation=1
.....

7FE0,0010,Pixel Data=131072

Patient Information

HIPPA Patient privacy

The Health Insurance Portability and Accountability Act of 1996 (HIPAA) protects the privacy and security of certain health information

<http://www.hhs.gov/hipaa/index.html>

0002,0000,File Meta Elements Group Len=148
0002,0001,File Meta Info Version=256
0002,0002,Media Storage SOP Class UID=1.2.840.10008.5.1.4.1.1.4.
0002,0003,Media Storage SOP Inst UID=0.0.0.0.
0002,0010,Transfer Syntax UID=1.2.840.10008.1.2.1.


...
0008,0060,Modality=MR
0008,0070,Manufacturer=GE MEDICAL SYSTEMS
0008,0080,Institution Name=1852796513
0008,0081,City Name=1852796513
0008,0090,Referring Physician's Name=1852796513
0008,0092,?=1852796513
0008,0201,?=-0500
0008,1010,Station Name=1852796513
0008,1030,Study Description=anon
0008,1035,Series Description=anon
0008,
0008,
0008,
0008,
0008,
0008,1000,Authorizing Diagnosis Description=1052796513
0009,1000,Manufacturer's Model Name=GENEPI SIGNA

Anonymized Patient Information

0010,0010,Patient's Name=anon
0010,0020,Patient ID=anon
0010,0030,Patient Date of Birth=00000000
0010,0032,Patient Birth Time=000000
0010,0040,Patient Sex=O
0010,1010,Patient Age=000Y

.....
0028,0010,Rows=256
0028,0011,Columns=256
0028,0030,Pixel Spacing=0.937500 0.937500
0028,0100,Bits Allocated=16
0028,0101,Bits Stored=16
0028,0102,High Bit=15
0028,0103,Pixel Representation=1
.....

7FE0,0010,Pixel Data=131072

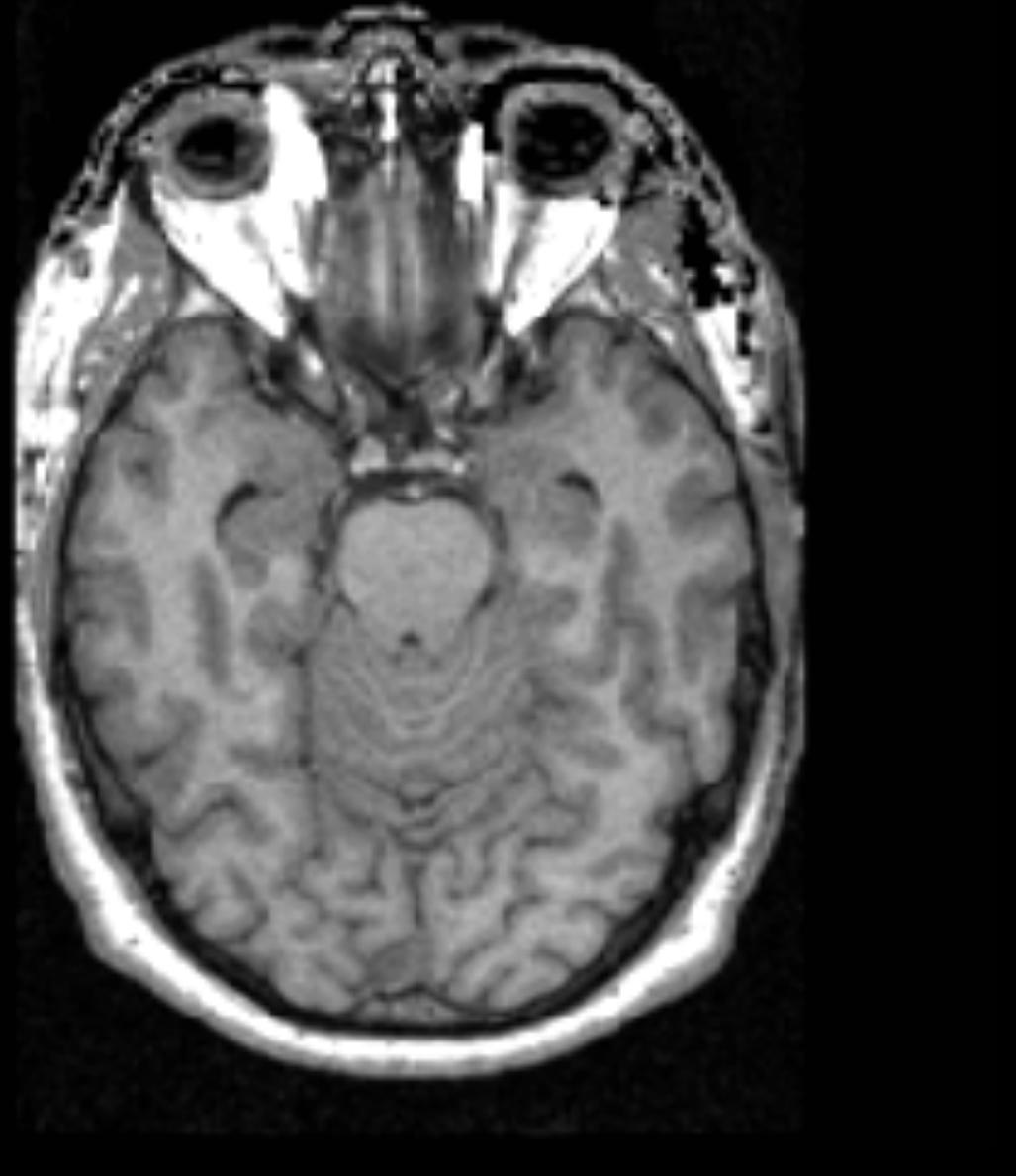

0002,0000,File Meta Elements Group Len=148
0002,0001,File Meta Info Version=256
0002,0002,Media Storage SOP Class UID=1.2.840.10008.5.1.4.1.1.4.
0002,0003,Media Storage SOP Inst UID=0.0.0.0.
0002,0010,Transfer Syntax UID=1.2.840.10008.1.2.1.


...

0008,0060,Modality=MR
0008,0070,Manufacturer=GE MEDICAL SYSTEMS
0008,0080,Institution Name=1852796513
0008,0081,City Name=1852796513
0008,0090,Referring Physician's Name=1852796513
0008,0092,?=1852796513
0008,0201,?=-0500
0008,1010,Station Name=1852796513
0008,1030,Study Description=anon
0008,103E,Series Description=anon
0008,1040,Institutional Dept. Name=1852796513
0008,1050,Performing Physician's Name=1852796513
0008,1060,Name Phys(s) Read Study=1852796513
0008,1070,Operator's Name=anon
0008,1080,Admitting Diagnosis Description=1852796513
0008,1090,Manufacturer's Model Name=GENESIS.SIGMA
0010,0100, Patient Age=0001
0010,0
0010,0
0010,0
0010,0
0010,1010, Patient Age=0001

Image Information

0028,0010,Rows=256
0028,0011,Columns=256
0028,0030,Pixel Spacing=0.937500 0.937500
0028,0100,Bits Allocated=16
0028,0101,Bits Stored=16
0028,0102,High Bit=15
0028,0103,Pixel Representation=1

7FE0,0010,Pixel Data=131072

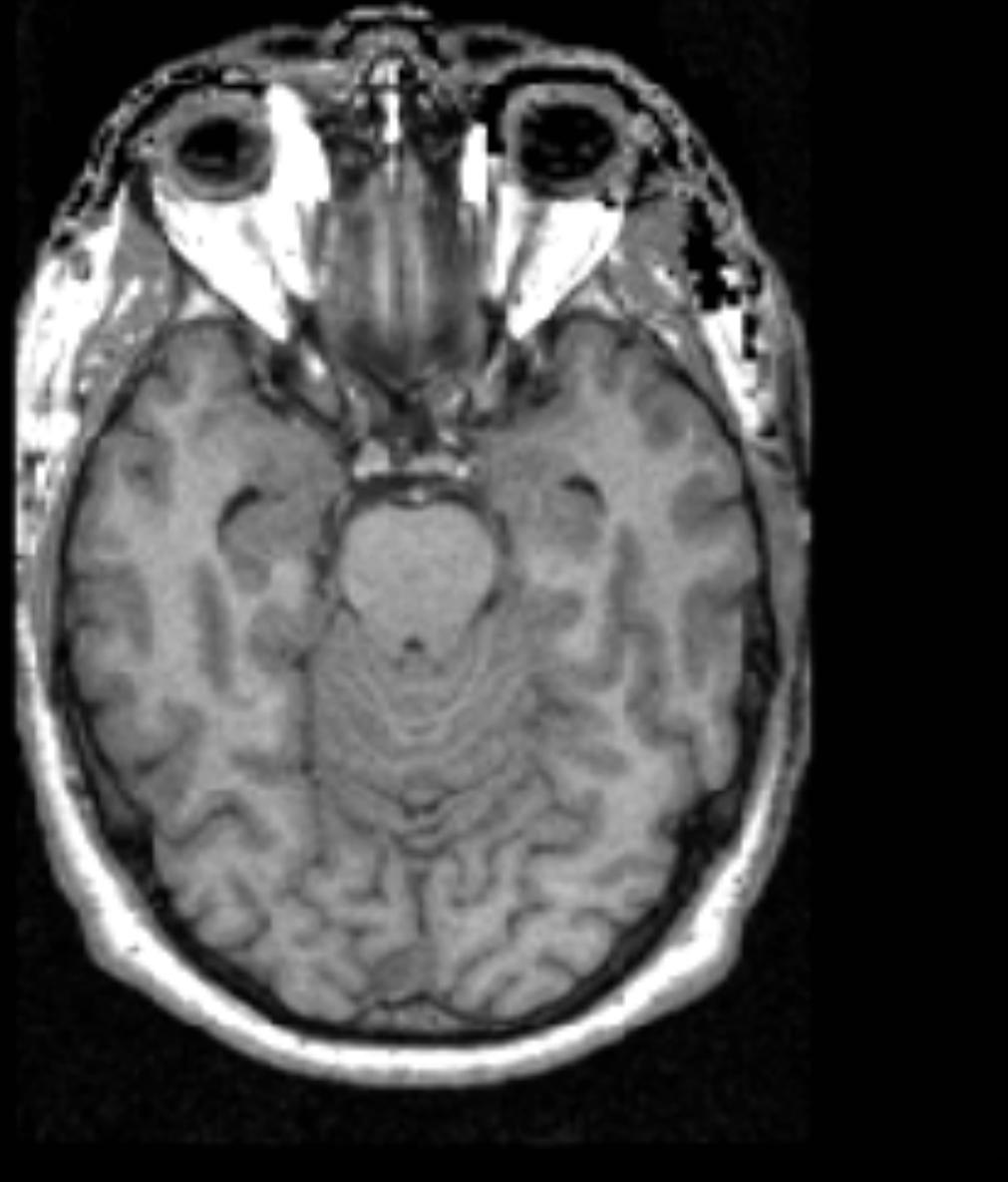

0002,0000,File Meta Elements Group Len=148
0002,0001,File Meta Info Version=256
0002,0002,Media Storage SOP Class UID=1.2.840.10008.5.1.4.1.1.4.
0002,0003,Media Storage SOP Inst UID=0.0.0.0.
0002,0010,Transfer Syntax UID=1.2.840.10008.1.2.1.

...

0008,0060,Modality=MR
0008,0070,Manufacturer=GE MEDICAL SYSTEMS
0008,0080,Institution Name=1852796513
0008,0081,City Name=1852796513
0008,0090,Referring Physician's Name=1852796513
0008,0092,?=1852796513
0008,0201,?=-0500
0008,1010,Station Name=1852796513
0008,1030,Study Description=anon
0008,103E,Series Description=anon
0008,1040,Institutional Dept. Name=1852796513
0008,1050,Performing Physician's Name=1852796513
0008,1060,Name Phys(s) Read Study=1852796513
0008,1070,Operator's Name=anon
0008,1080,Admitting Diagnosis Description=1852796513
0008,1090,Manufacturer's Model Name=GENESIS.SIGMA


0010,0010,Patient's Name=anon
0010,0020,Patient ID=anon
0010,0030,Patient Date of Birth=00000000
0010,0032,Patient Birth Time=000000
0010,0040,Patient Sex=O
0010,1010,Patient Age=000Y

.....
0028,0010,Rows=256
0028,0011,Cols=256
0028,0012,Pixel Depth=8
0028,0013,Pixel Representation=1
0028,0014,Pixel Type=1
0028,0015,Pixel Spacing=1.000000
0028,0016,Pixel Pitch=1.000000
0028,0017,Pixel Aspect Ratio=1.000000
0028,0018,Image Resolution=1.000000
0028,0019,Image Resolution Units=MILLIMETER
0028,001A,Image Resolution Dimensions=2
0028,001B,Image Resolution Interpolate=0
0028,001C,Image Resolution Interpolate X=0
0028,001D,Image Resolution Interpolate Y=0
0028,001E,Image Resolution Interpolate Z=0
0028,001F,Image Resolution Interpolate X Units=MILLIMETER
0028,0020,Image Resolution Interpolate Y Units=MILLIMETER
0028,0021,Image Resolution Interpolate Z Units=MILLIMETER
0028,0022,Image Resolution Interpolate X Min=0
0028,0023,Image Resolution Interpolate Y Min=0
0028,0024,Image Resolution Interpolate Z Min=0
0028,0025,Image Resolution Interpolate X Max=1
0028,0026,Image Resolution Interpolate Y Max=1
0028,0027,Image Resolution Interpolate Z Max=1
0028,0028,Image Resolution Interpolate X Step=1
0028,0029,Image Resolution Interpolate Y Step=1
0028,002A,Image Resolution Interpolate Z Step=1
0028,002B,Image Resolution Interpolate X Min Value=0
0028,002C,Image Resolution Interpolate Y Min Value=0
0028,002D,Image Resolution Interpolate Z Min Value=0
0028,002E,Image Resolution Interpolate X Max Value=1
0028,002F,Image Resolution Interpolate Y Max Value=1
0028,0030,Image Resolution Interpolate Z Max Value=1
0028,0031,Image Resolution Interpolate X Step Value=1
0028,0032,Image Resolution Interpolate Y Step Value=1
0028,0033,Image Resolution Interpolate Z Step Value=1
0028,0034,Image Resolution Interpolate X Min Value=0
0028,0035,Image Resolution Interpolate Y Min Value=0
0028,0036,Image Resolution Interpolate Z Min Value=0
0028,0037,Image Resolution Interpolate X Max Value=1
0028,0038,Image Resolution Interpolate Y Max Value=1
0028,0039,Image Resolution Interpolate Z Max Value=1
0028,003A,Image Resolution Interpolate X Step Value=1
0028,003B,Image Resolution Interpolate Y Step Value=1
0028,003C,Image Resolution Interpolate Z Step Value=1
0028,003D,Image Resolution Interpolate X Min Value=0
0028,003E,Image Resolution Interpolate Y Min Value=0
0028,003F,Image Resolution Interpolate Z Min Value=0
0028,0040,Image Resolution Interpolate X Max Value=1
0028,0041,Image Resolution Interpolate Y Max Value=1
0028,0042,Image Resolution Interpolate Z Max Value=1
0028,0043,Image Resolution Interpolate X Step Value=1
0028,0044,Image Resolution Interpolate Y Step Value=1
0028,0045,Image Resolution Interpolate Z Step Value=1
0028,0046,Image Resolution Interpolate X Min Value=0
0028,0047,Image Resolution Interpolate Y Min Value=0
0028,0048,Image Resolution Interpolate Z Min Value=0
0028,0049,Image Resolution Interpolate X Max Value=1
0028,004A,Image Resolution Interpolate Y Max Value=1
0028,004B,Image Resolution Interpolate Z Max Value=1
0028,004C,Image Resolution Interpolate X Step Value=1
0028,004D,Image Resolution Interpolate Y Step Value=1
0028,004E,Image Resolution Interpolate Z Step Value=1
0028,004F,Image Resolution Interpolate X Min Value=0
0028,0050,Image Resolution Interpolate Y Min Value=0
0028,0051,Image Resolution Interpolate Z Min Value=0
0028,0052,Image Resolution Interpolate X Max Value=1
0028,0053,Image Resolution Interpolate Y Max Value=1
0028,0054,Image Resolution Interpolate Z Max Value=1
0028,0055,Image Resolution Interpolate X Step Value=1
0028,0056,Image Resolution Interpolate Y Step Value=1
0028,0057,Image Resolution Interpolate Z Step Value=1
0028,0058,Image Resolution Interpolate X Min Value=0
0028,0059,Image Resolution Interpolate Y Min Value=0
0028,005A,Image Resolution Interpolate Z Min Value=0
0028,005B,Image Resolution Interpolate X Max Value=1
0028,005C,Image Resolution Interpolate Y Max Value=1
0028,005D,Image Resolution Interpolate Z Max Value=1
0028,005E,Image Resolution Interpolate X Step Value=1
0028,005F,Image Resolution Interpolate Y Step Value=1
0028,0060,Image Resolution Interpolate Z Step Value=1
0028,0061,Image Resolution Interpolate X Min Value=0
0028,0062,Image Resolution Interpolate Y Min Value=0
0028,0063,Image Resolution Interpolate Z Min Value=0
0028,0064,Image Resolution Interpolate X Max Value=1
0028,0065,Image Resolution Interpolate Y Max Value=1
0028,0066,Image Resolution Interpolate Z Max Value=1
0028,0067,Image Resolution Interpolate X Step Value=1
0028,0068,Image Resolution Interpolate Y Step Value=1
0028,0069,Image Resolution Interpolate Z Step Value=1
0028,006A,Image Resolution Interpolate X Min Value=0
0028,006B,Image Resolution Interpolate Y Min Value=0
0028,006C,Image Resolution Interpolate Z Min Value=0
0028,006D,Image Resolution Interpolate X Max Value=1
0028,006E,Image Resolution Interpolate Y Max Value=1
0028,006F,Image Resolution Interpolate Z Max Value=1
0028,0070,Image Resolution Interpolate X Step Value=1
0028,0071,Image Resolution Interpolate Y Step Value=1
0028,0072,Image Resolution Interpolate Z Step Value=1
0028,0073,Image Resolution Interpolate X Min Value=0
0028,0074,Image Resolution Interpolate Y Min Value=0
0028,0075,Image Resolution Interpolate Z Min Value=0
0028,0076,Image Resolution Interpolate X Max Value=1
0028,0077,Image Resolution Interpolate Y Max Value=1
0028,0078,Image Resolution Interpolate Z Max Value=1
0028,0079,Image Resolution Interpolate X Step Value=1
0028,0080,Image Resolution Interpolate Y Step Value=1
0028,0081,Image Resolution Interpolate Z Step Value=1
0028,0082,Image Resolution Interpolate X Min Value=0
0028,0083,Image Resolution Interpolate Y Min Value=0
0028,0084,Image Resolution Interpolate Z Min Value=0
0028,0085,Image Resolution Interpolate X Max Value=1
0028,0086,Image Resolution Interpolate Y Max Value=1
0028,0087,Image Resolution Interpolate Z Max Value=1
0028,0088,Image Resolution Interpolate X Step Value=1
0028,0089,Image Resolution Interpolate Y Step Value=1
0028,0090,Image Resolution Interpolate Z Step Value=1
0028,0091,Image Resolution Interpolate X Min Value=0
0028,0092,Image Resolution Interpolate Y Min Value=0
0028,0093,Image Resolution Interpolate Z Min Value=0
0028,0094,Image Resolution Interpolate X Max Value=1
0028,0095,Image Resolution Interpolate Y Max Value=1
0028,0096,Image Resolution Interpolate Z Max Value=1
0028,0097,Image Resolution Interpolate X Step Value=1
0028,0098,Image Resolution Interpolate Y Step Value=1
0028,0099,Image Resolution Interpolate Z Step Value=1
0028,0010,Pixel Data=131072

Pixel Data

Standard DICOM Instances

- **DICOM images** produced by imaging equipment: single DICOM MRI slices, DICOM CT volume, DICOM multi-frame objects (fMRI experiment)
- **DICOM Segmentation Object (SEG)**: voxels labelled in regions of interest (ROIs)
- **DICOM Structured Report (SR)**: clinical information (e.g. diagnosis, pathology, surgery, etc.), measurements computed from segmented ROIs

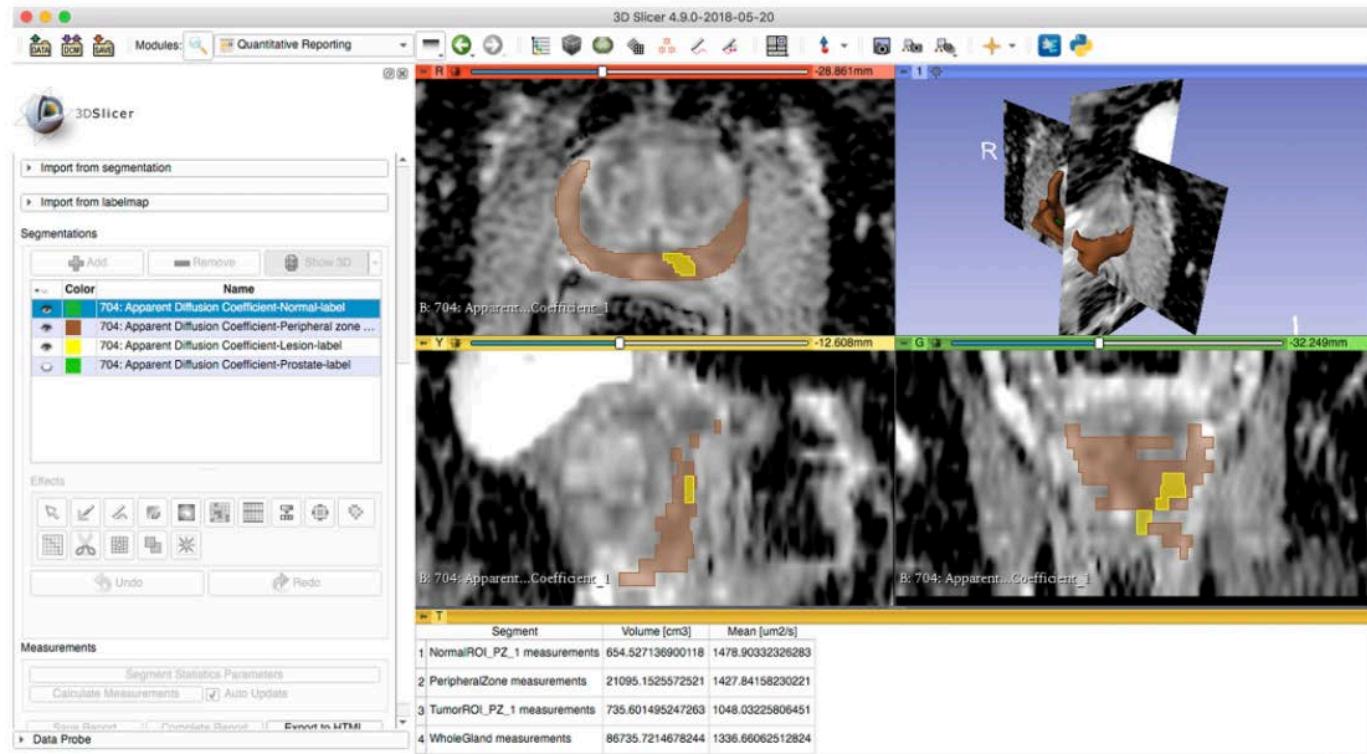

Examples of Standard DICOM Instances

DICOM Images

Prostate MRI Images

DICOM Segmented Structures

Prostate gland, peripheral zone, lesion, normal tissue

Segmented structure	SegmentedPropertyCategoryCodeSequence	SegmentedPropertyTypeCodeSequence	AnatomicRegionSequence
Prostate gland	("T-D000A", "SRT", "Anatomical Structure")	("T-9200B", "SRT", "Prostate")	NA
Peripheral zone of the prostate gland	("T-D000A", "SRT", "Anatomical Structure")	("T-D05E4", "SRT", "Peripheral zone of the prostate")	NA
Lesion in the peripheral zone of the prostate gland	("M-01000", "SRT", "Morphologically Altered Structure")	("M-01100", "SRT", "Lesion")	("T-D05E4", "SRT", "Peripheral zone of the prostate")
Normal tissue in the peripheral zone of the prostate gland	("T-D0050", "SRT", "Tissue")	("G-A460", "SRT", "Normal")	("T-D05E4", "SRT", "Peripheral zone of the prostate")



Measured structure	Finding	Finding Site
Prostate gland *	(T-F6078, SRT, "Entire Gland")	("T-9200B", "SRT", "Prostate")
Peripheral zone of the prostate gland	(R-404A4, SRT, "Entire")	("T-D05E4", "SRT", "Peripheral zone of the prostate")
Lesion in the peripheral zone of the prostate gland	(R-42037, SRT, "Abnormal")	("T-D05E4", "SRT", "Peripheral zone of the prostate")
Normal tissue of peripheral zone of the prostate gland	("G-A460", "SRT", "Normal")	("T-D05E4", "SRT", "Peripheral zone of the prostate")

Measurements

Volumes of prostate gland, peripheral zone, lesion, normal tissue

DICOM Terminology

- DICOM provides support for annotation and segmentation objects
- DICOM Structured Reports enable provenance tracking
- DICOM derived data can be stored on a DICOM server or on other archive (e.g. TCIA) with permission and is compliant with FAIR principles

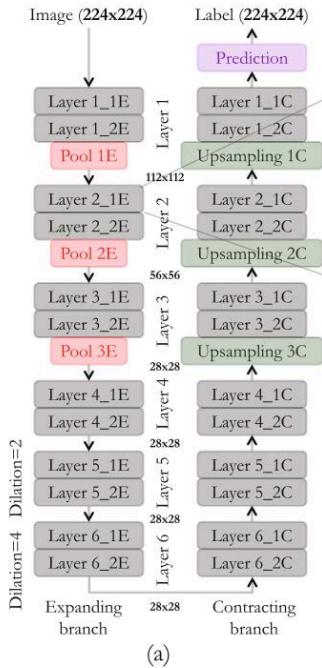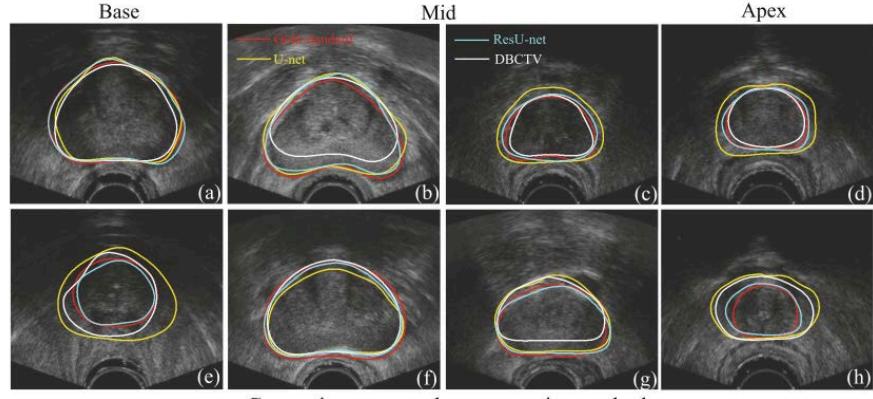


DICOM for AI Studies

- DICOM defines **syntax rules** and **vocabularies** that enable easy extraction of knowledge from the data
- The DICOM framework for medical imaging data management enables the automation of cohort formation and maximizes the interoperability of the data for AI studies



Part 2: Loading and Visualizing DICOM data in 3D Slicer version 5.0

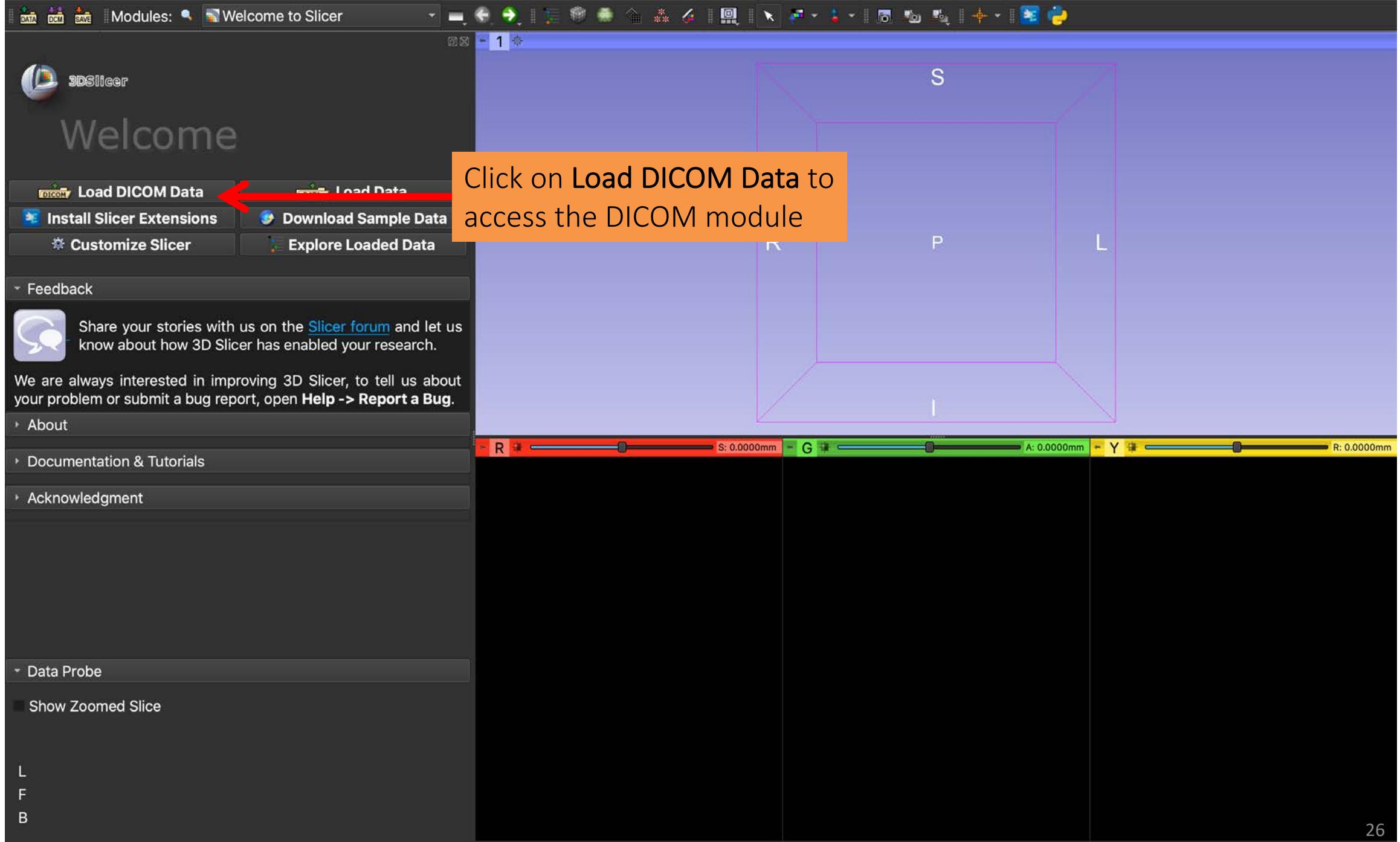

The DICOM database follows the DICOM model hierarchy organized in patients, studies and series

Patients

Studies

Series

DICOM Data Reader Warnings

Uncheck All Examine Load

DATA DCM SAVE Modules: DICOM DICOM database

3DSlicer Help & Acknowledgement

Import DICOM files Show DICOM database

Loaded data Node

DICOM networking DICOM database settings

Advanced 27


Sonia Pujol PhD, All Rights Reserved



Dataset #1

Torso CT

Loading a DICOM dataset in Slicer

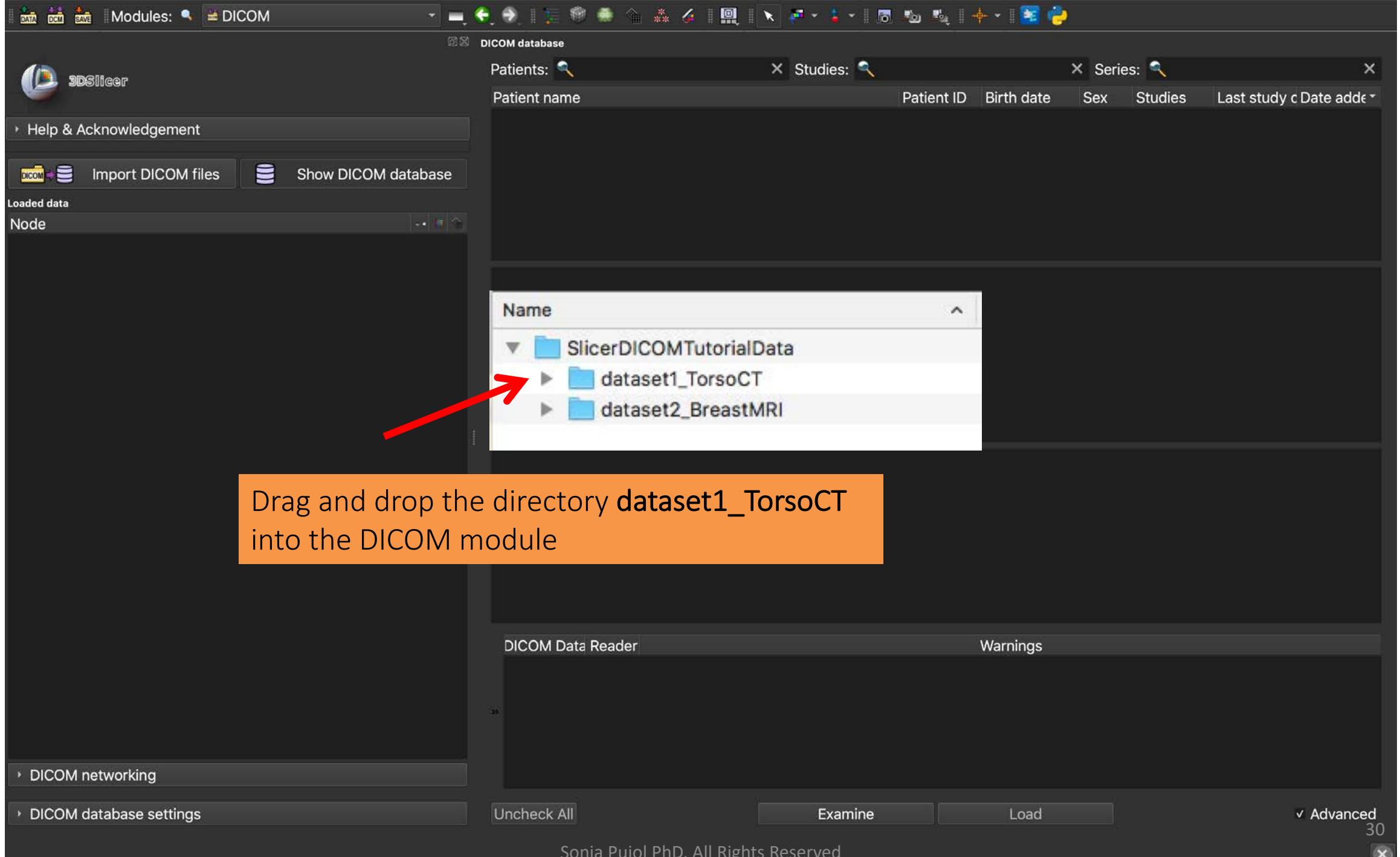

Click on patient1 in the list of patients

Slicer displays the corresponding study and series

The screenshot shows the 3DSlicer interface with the DICOM module open. In the top right, the 'DICOM database' panel is visible, displaying a table of patients, studies, and series. A red arrow points from the text 'Click on patient1 in the list of patients' to the row for 'patient1'. Another red arrow points from the text 'Slicer displays the corresponding study and series' to the 'Studies' and 'Series' tables below. The 'Studies' table shows one study for 'patient1_ID' with details: Study date 20050601, Study ID 6936864, and Study description CT Thorax Abdomen. The 'Series' table shows one series for the study with details: Series # 6, Series description CT_Thorax_Abdomen, Modality CT, Size 512x512, Count 291, and Date added 202...843. The bottom right corner of the slide indicates page 31.

Patient name	Patient ID	Birth date	Sex	Studies	Last study	Date add
patient1	patient1_ID			1	Wed Jun 1 2005	2020...842

Study date	Study ID	Study description	Series	Date added
20050601	6936864	CT Thorax Abdomen	1	202...843

Series #	Series description	Modality	Size	Count	Date added
6	CT_Thorax_Abdomen	CT	512x512	291	202...843

DICOM Data Reader Warnings
6: ... Scalar Volume

DICOM networking DICOM database settings Uncheck All Examine Load Advanced 31

3DSlicer DICOM Modules: DICOM

DICOM database

Patients: Studies: Series:

Patient name	Patient ID	Birth date	Sex	Studies	Last study	Date added
patient1	patient1_ID			1	Wed Jun 1 2020...	2005 842

Import DICOM files Show DICOM database

Loaded data Node

Study date Study ID Study description Series Date added

20050601	6936864	CT Thorax Abdomen	1	202...843
----------	---------	-------------------	---	-----------

Series # Series description Modality Size Count Date added

6	CT_Thorax_Abdomen	CT	512x512	291	202...843
---	-------------------	----	---------	-----	-----------

DICOM Data Reader Warnings

6: ... Scalar Volume

Click on Examine

Uncheck All Examine Load Advanced 32

DICOM networking DICOM database settings

Sonia Pujol PhD, All Rights Reserved

3DSlicer DICOM Modules: DICOM

DICOM database

Patients: Studies: Series:

Patient name	Patient ID	Birth date	Sex	Studies	Last study	Date added
patient1	patient1_ID			1	Wed Jun 1 2005	2020...842

Import DICOM files Show DICOM database

Loaded data Node

Study date Study ID Study description Series Date added

20050601	6936864	CT Thorax Abdomen	1	202...843
----------	---------	-------------------	---	-----------

Series # Series description Modality Size Count Date added

6	CT_Thorax_Abdomen	CT	512x512	291	202...843
---	-------------------	----	---------	-----	-----------

DICOM Data Reader Warnings


6: ... Scalar Volume

Click on Load


Uncheck All Examine Load Advanced

DICOM networking DICOM database settings


Sonia Pujol PhD, All Rights Reserved

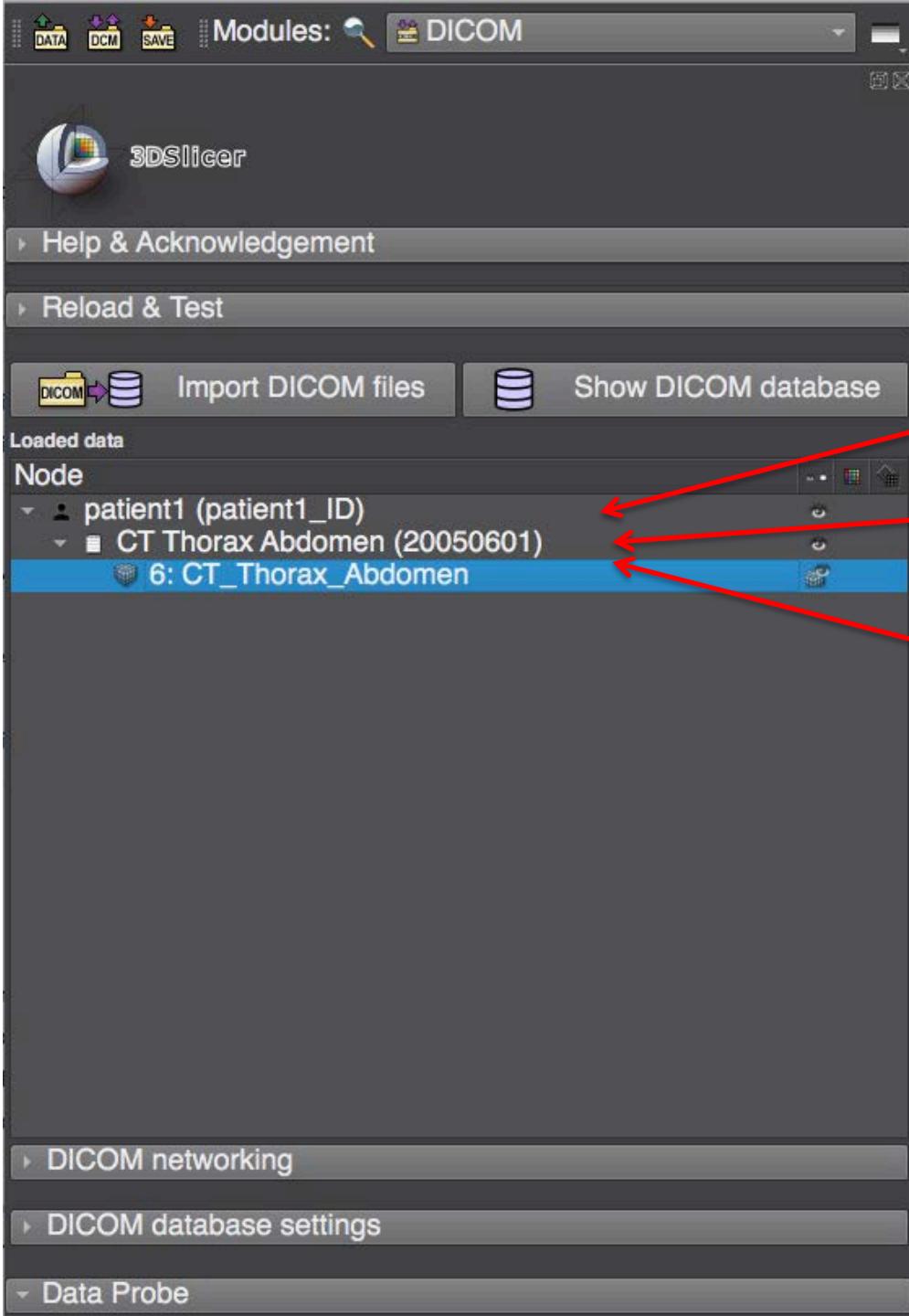


The DICOM dataset is loaded into Slicer as a patient-study-series DICOM hierarchy

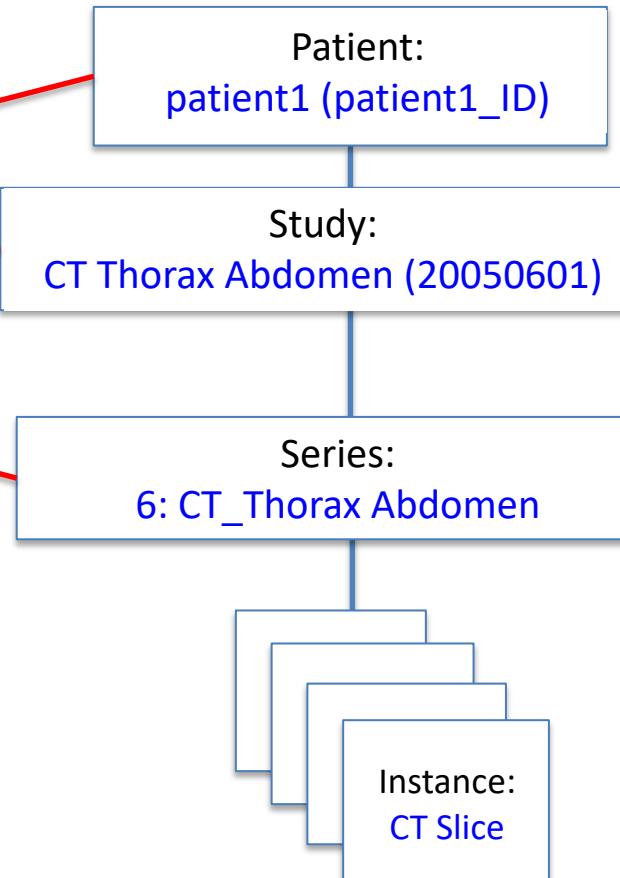



Visualizing a DICOM dataset in Slicer

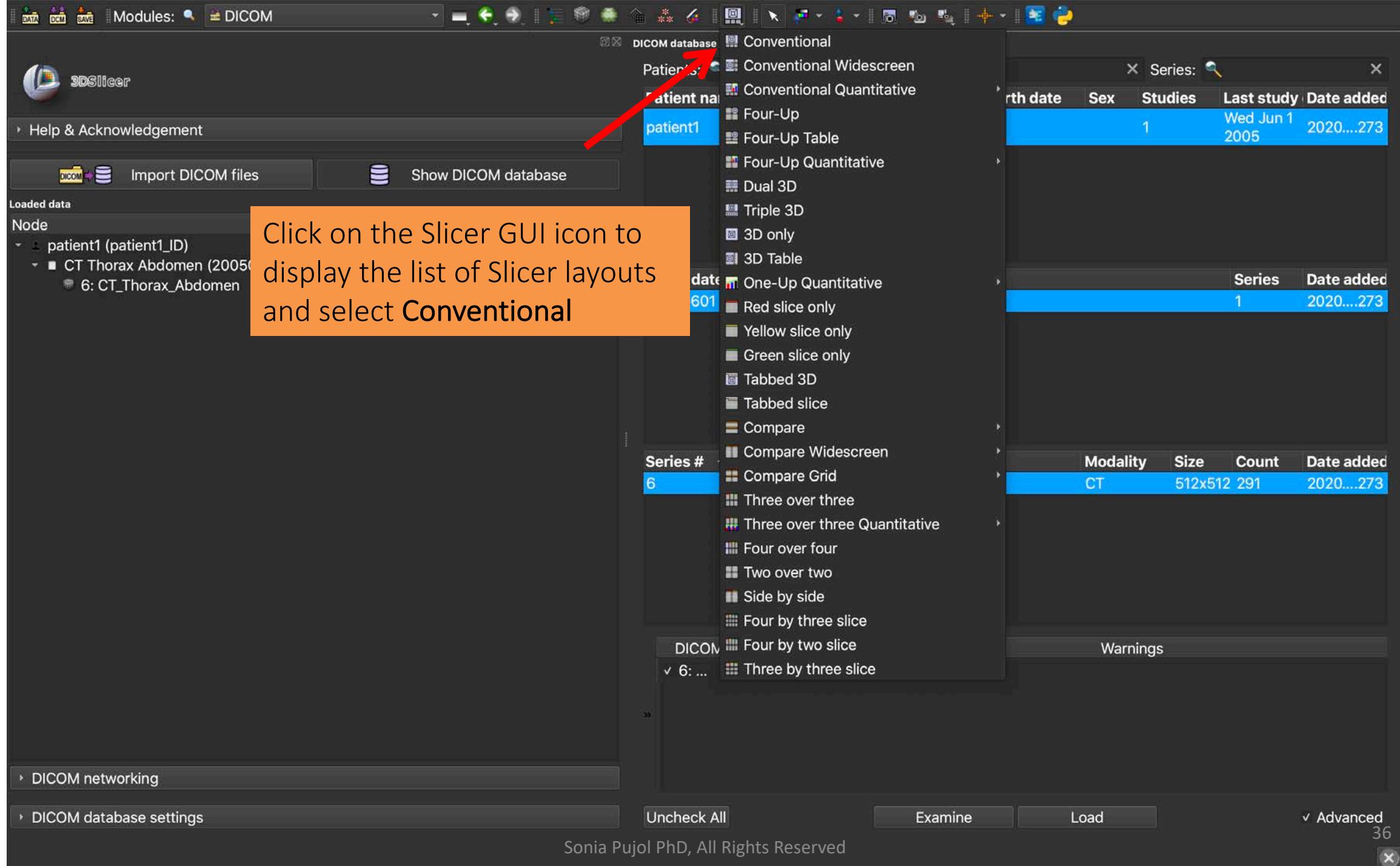



Click on the Slicer GUI icon to display the list of Slicer layouts and select Conventional

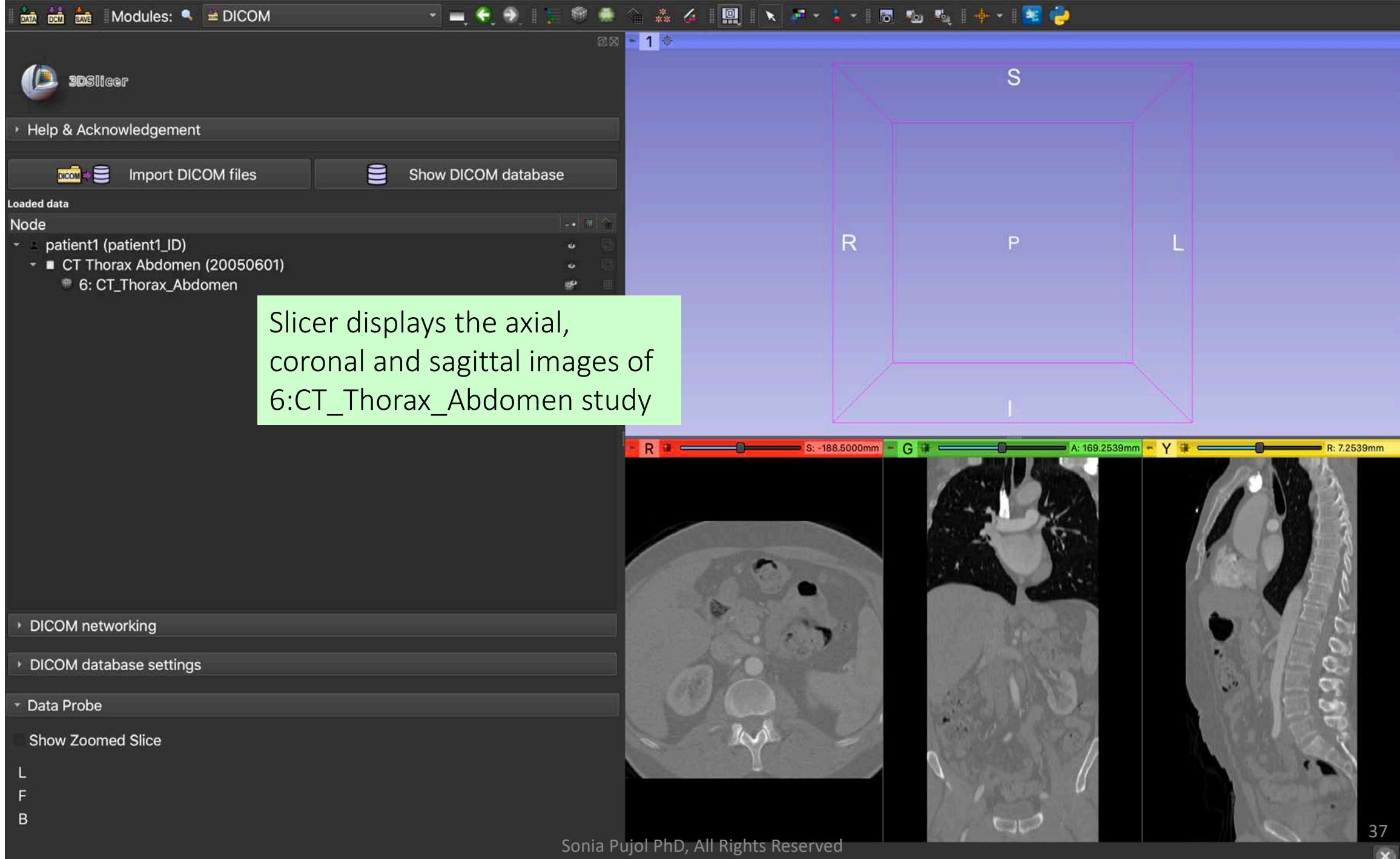


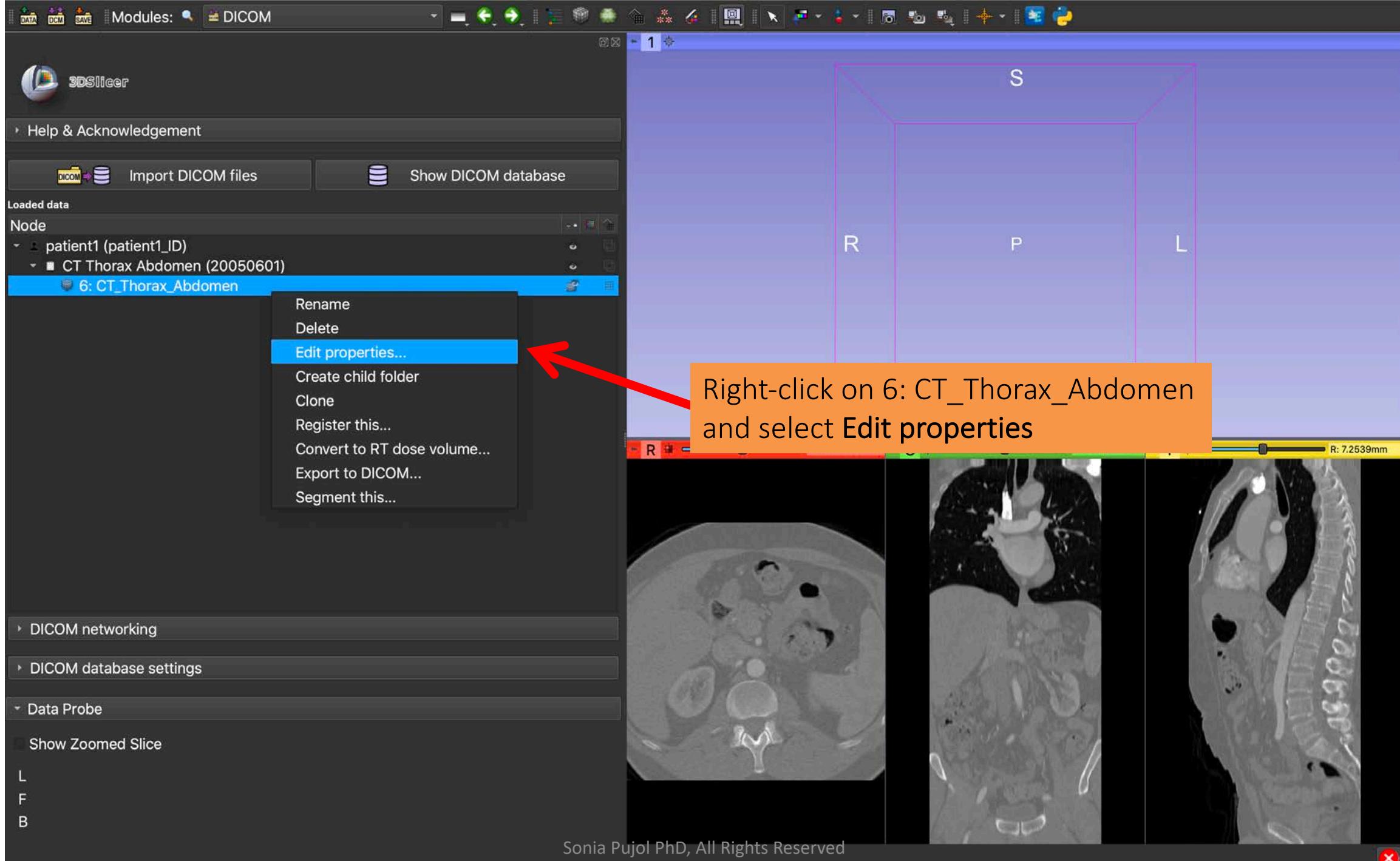


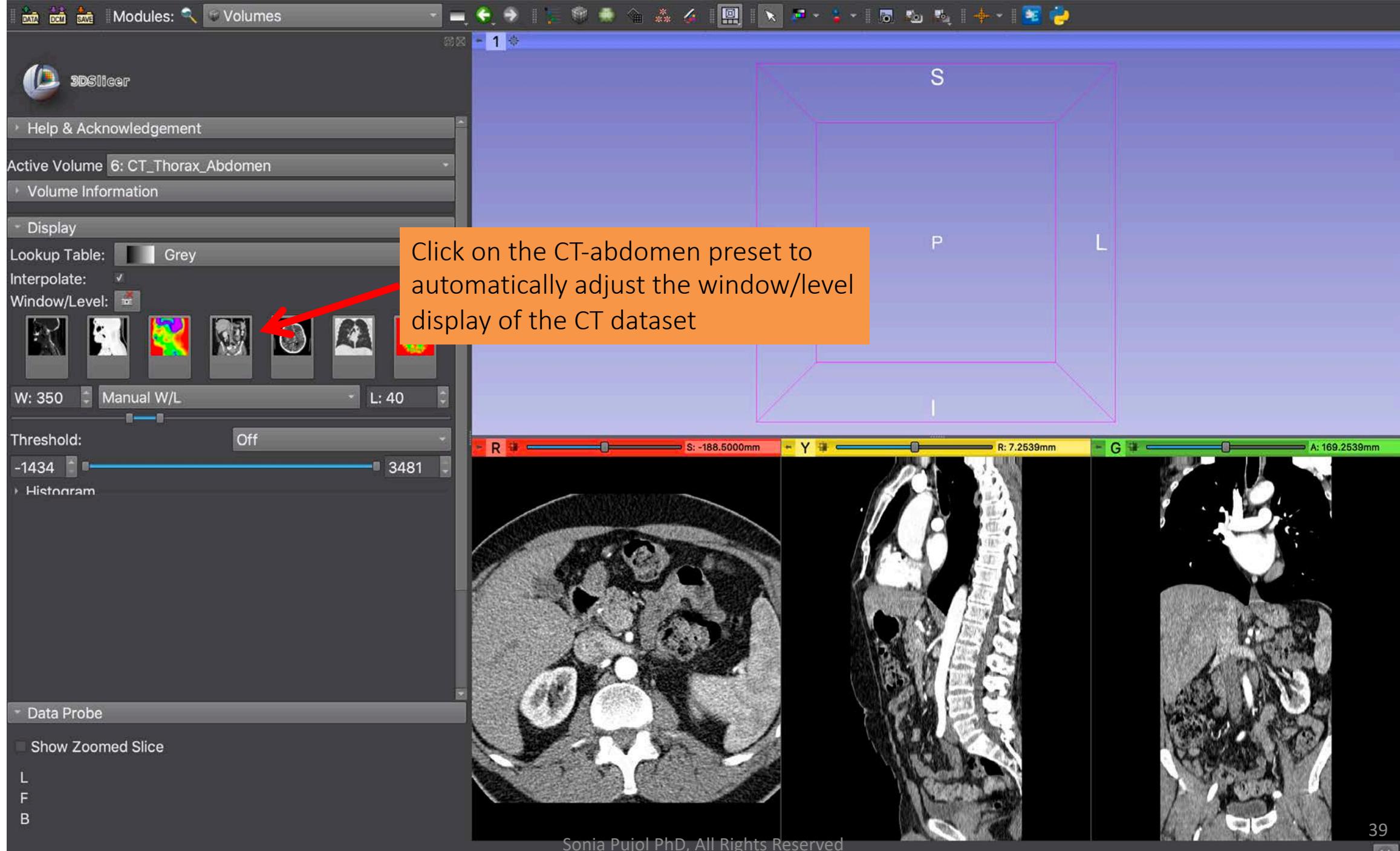


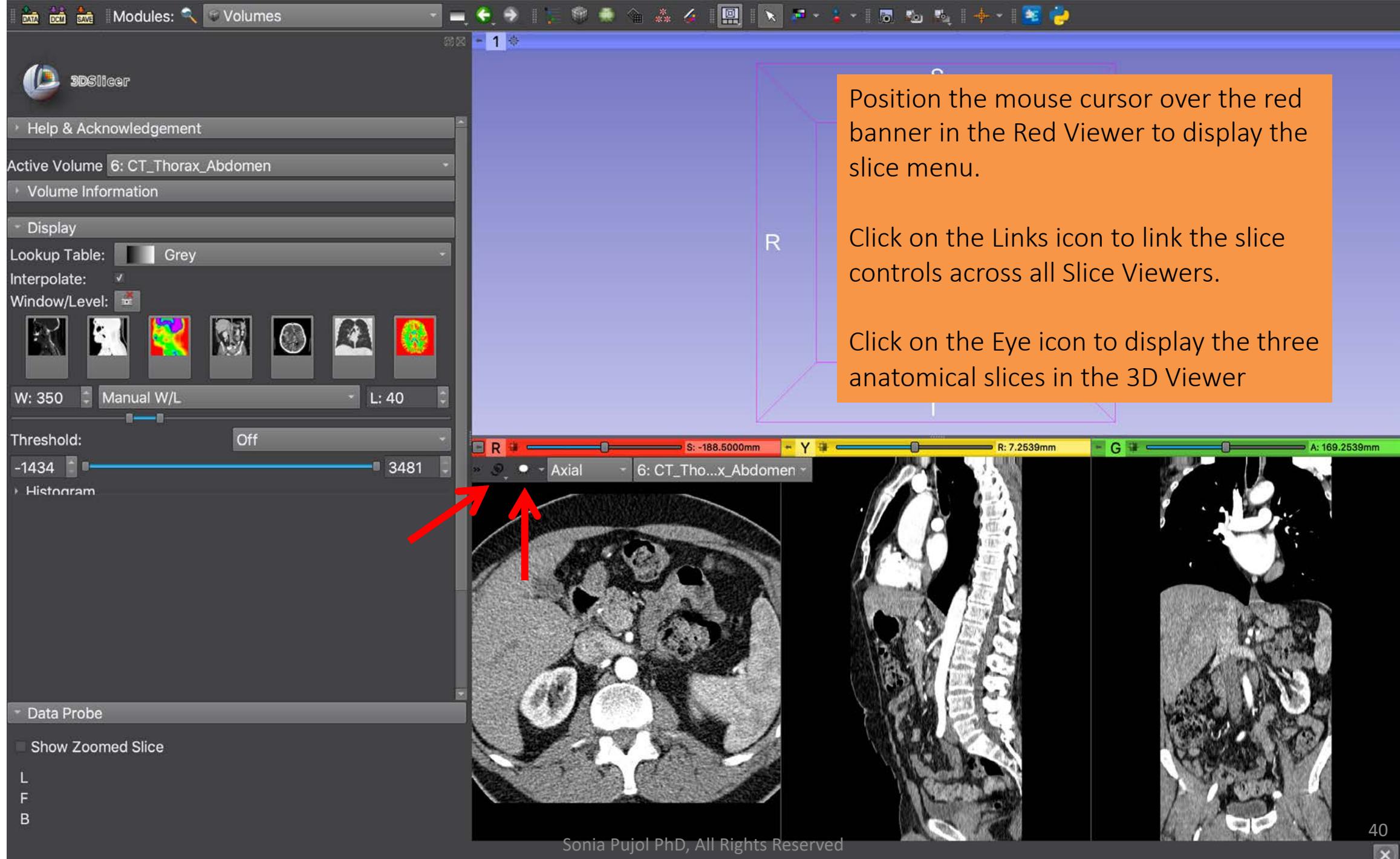


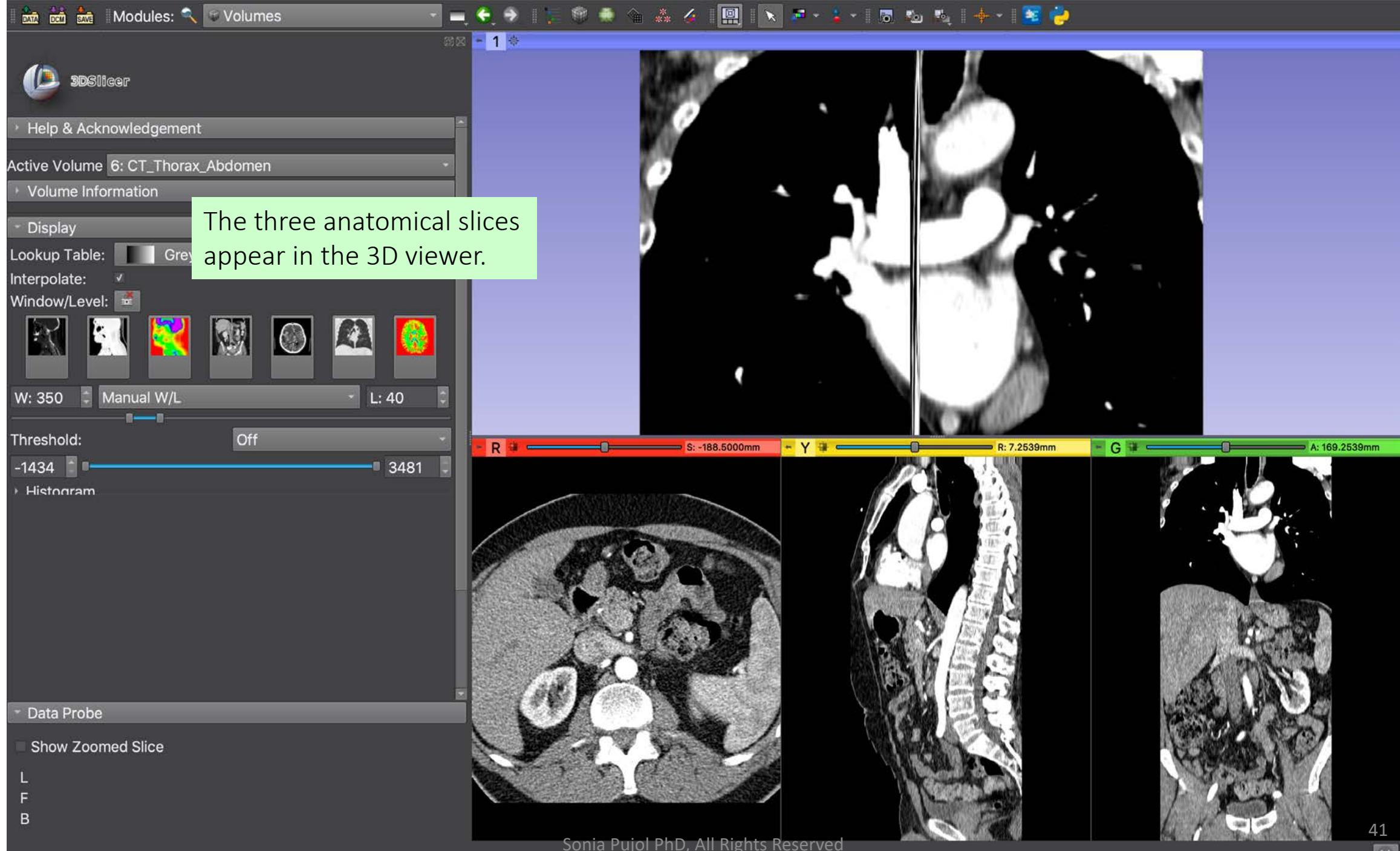


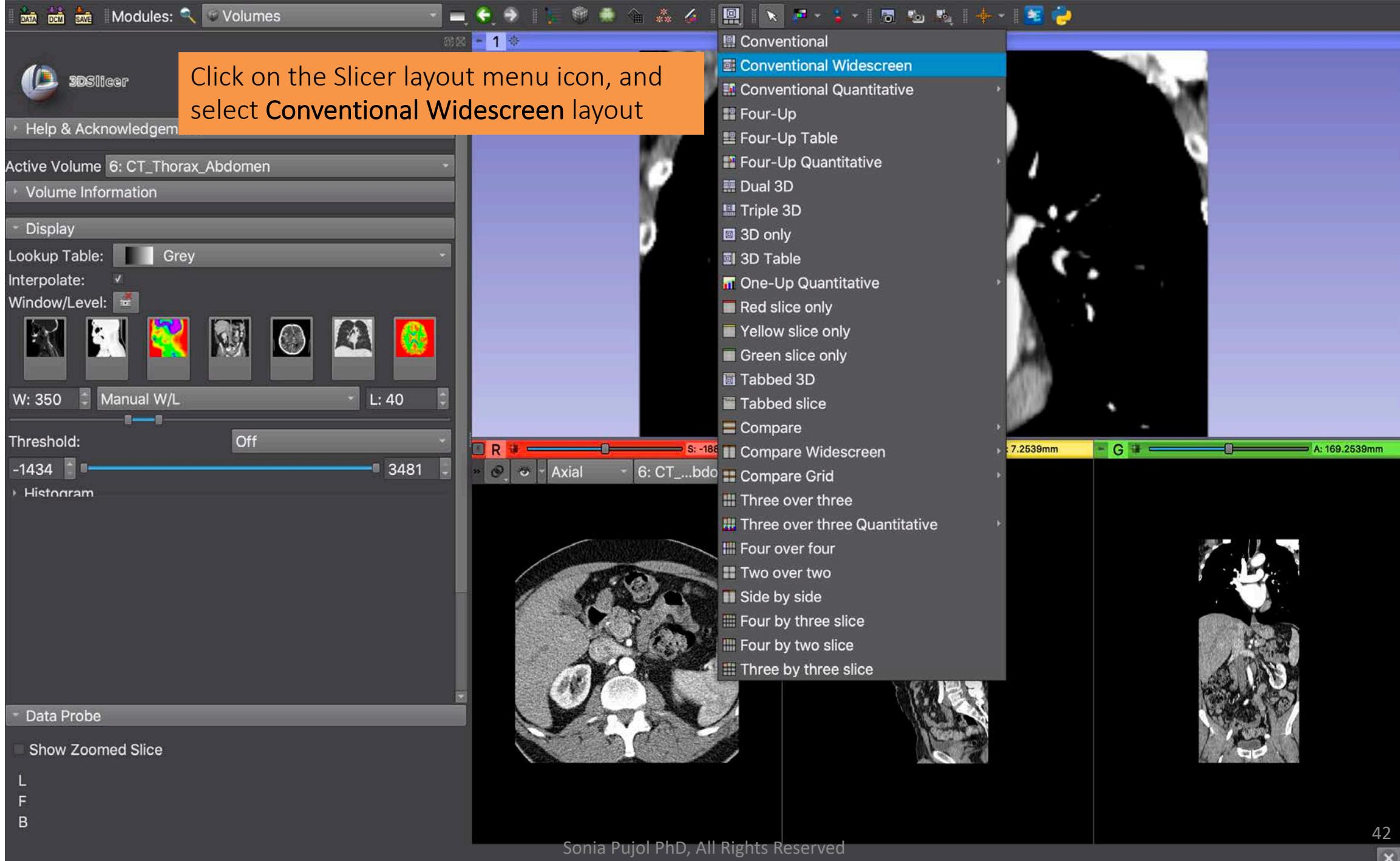


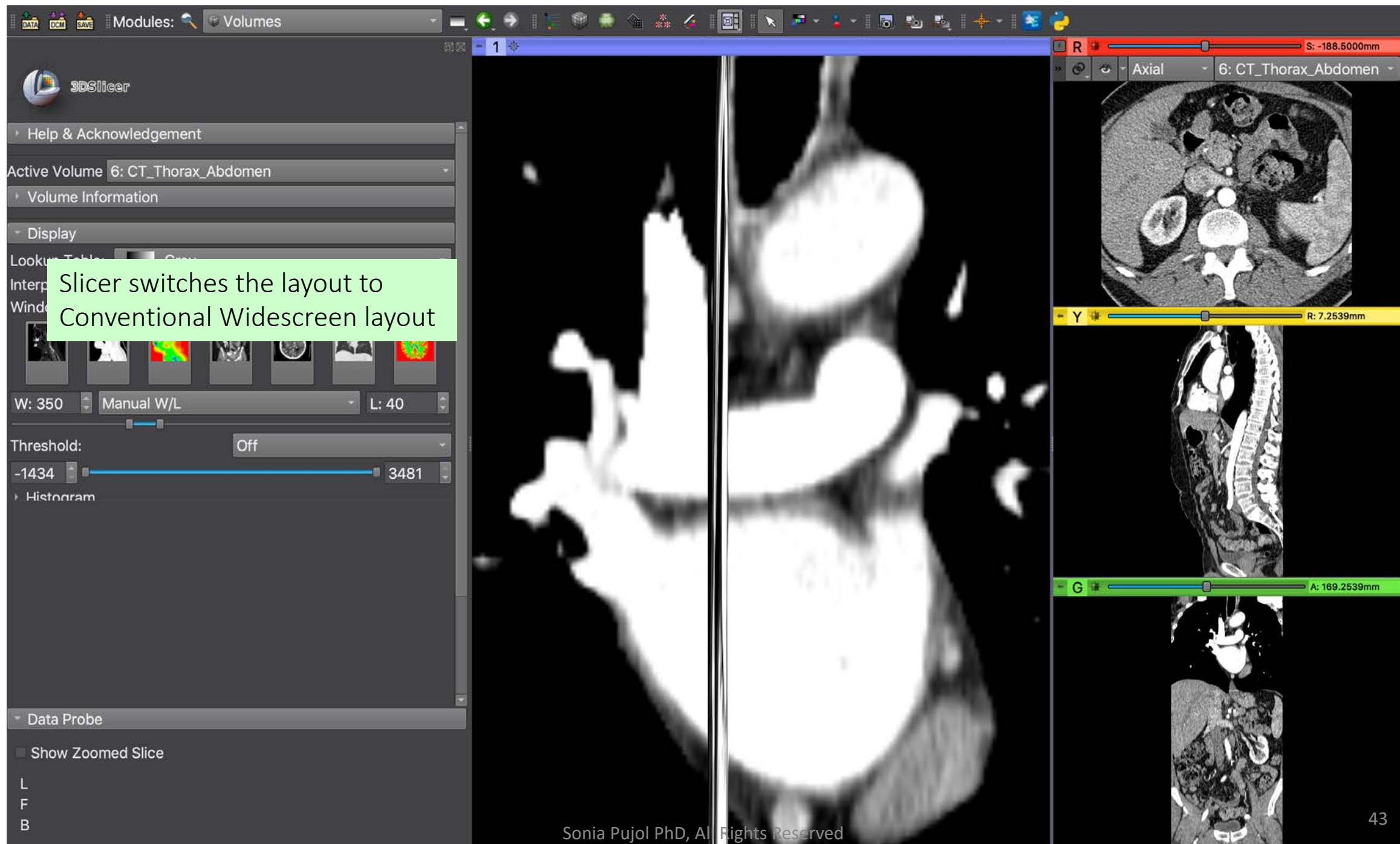


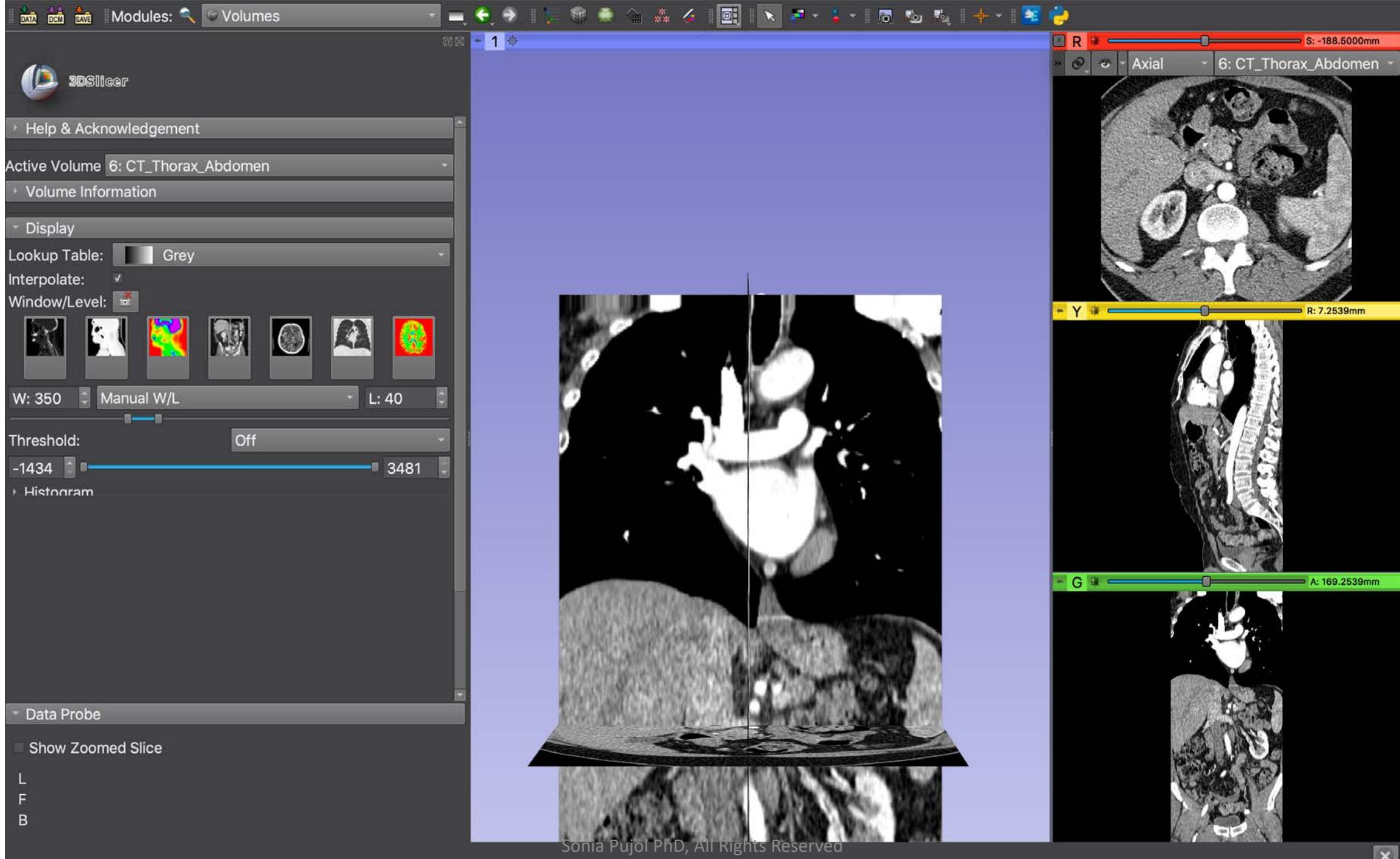

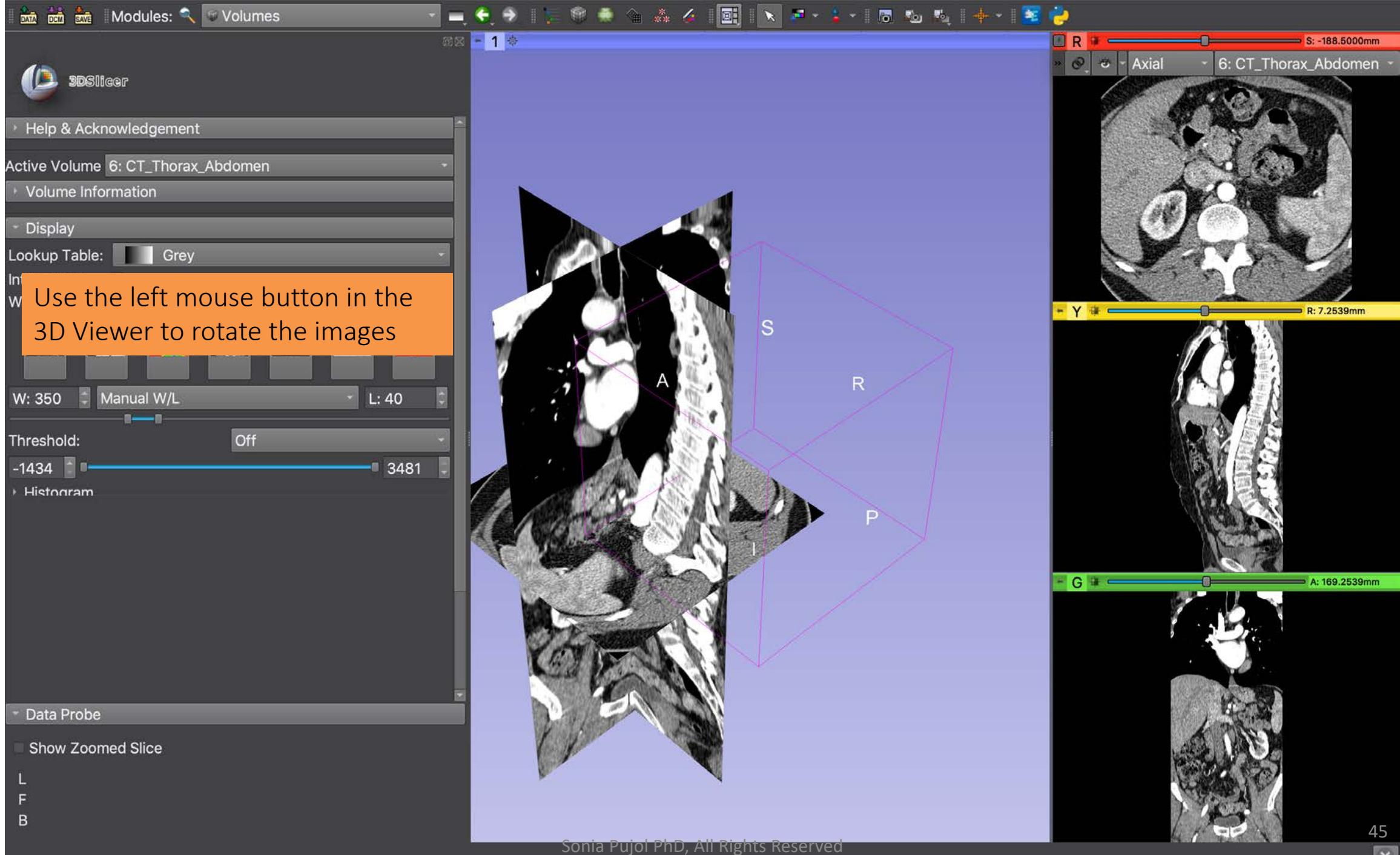

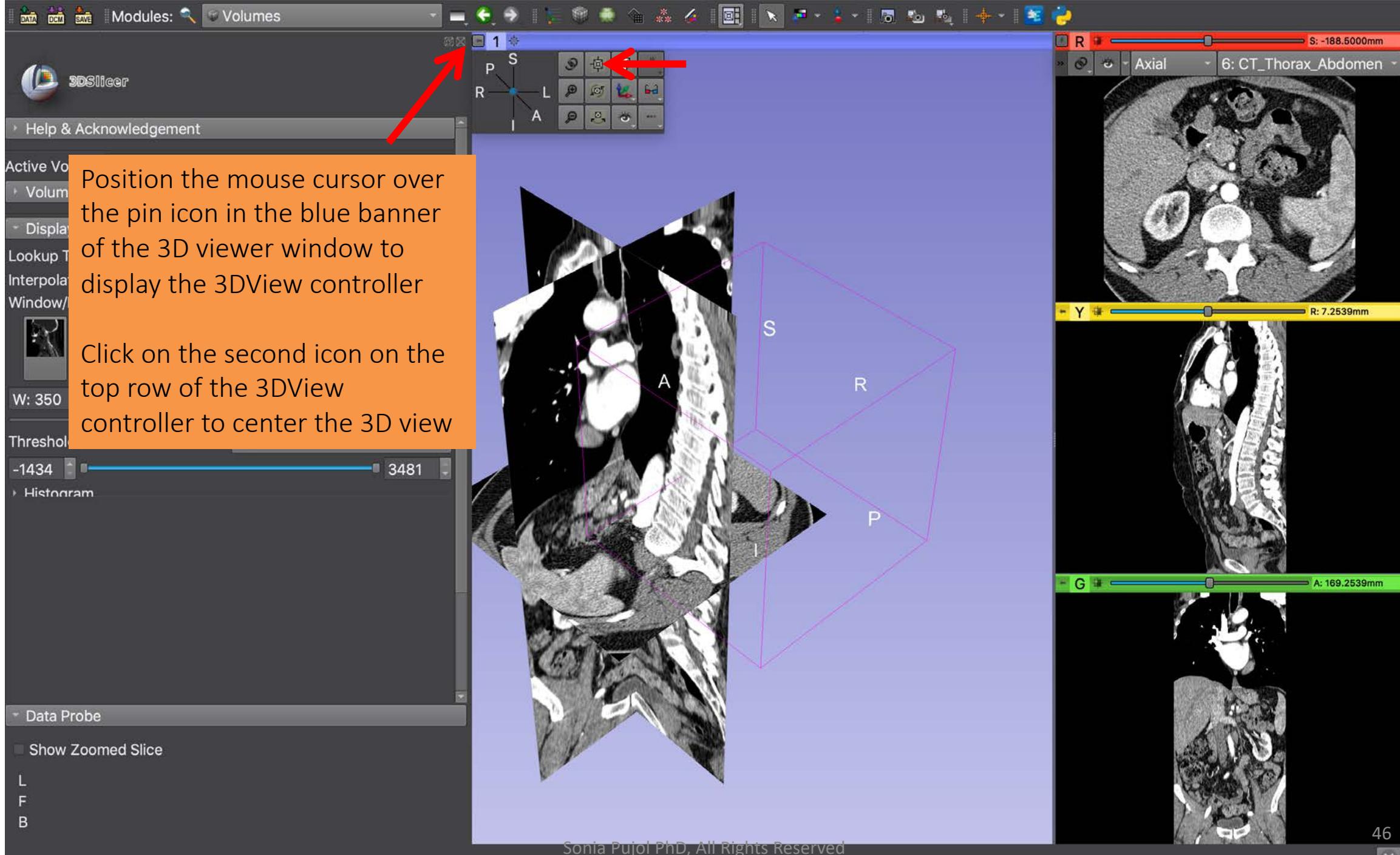

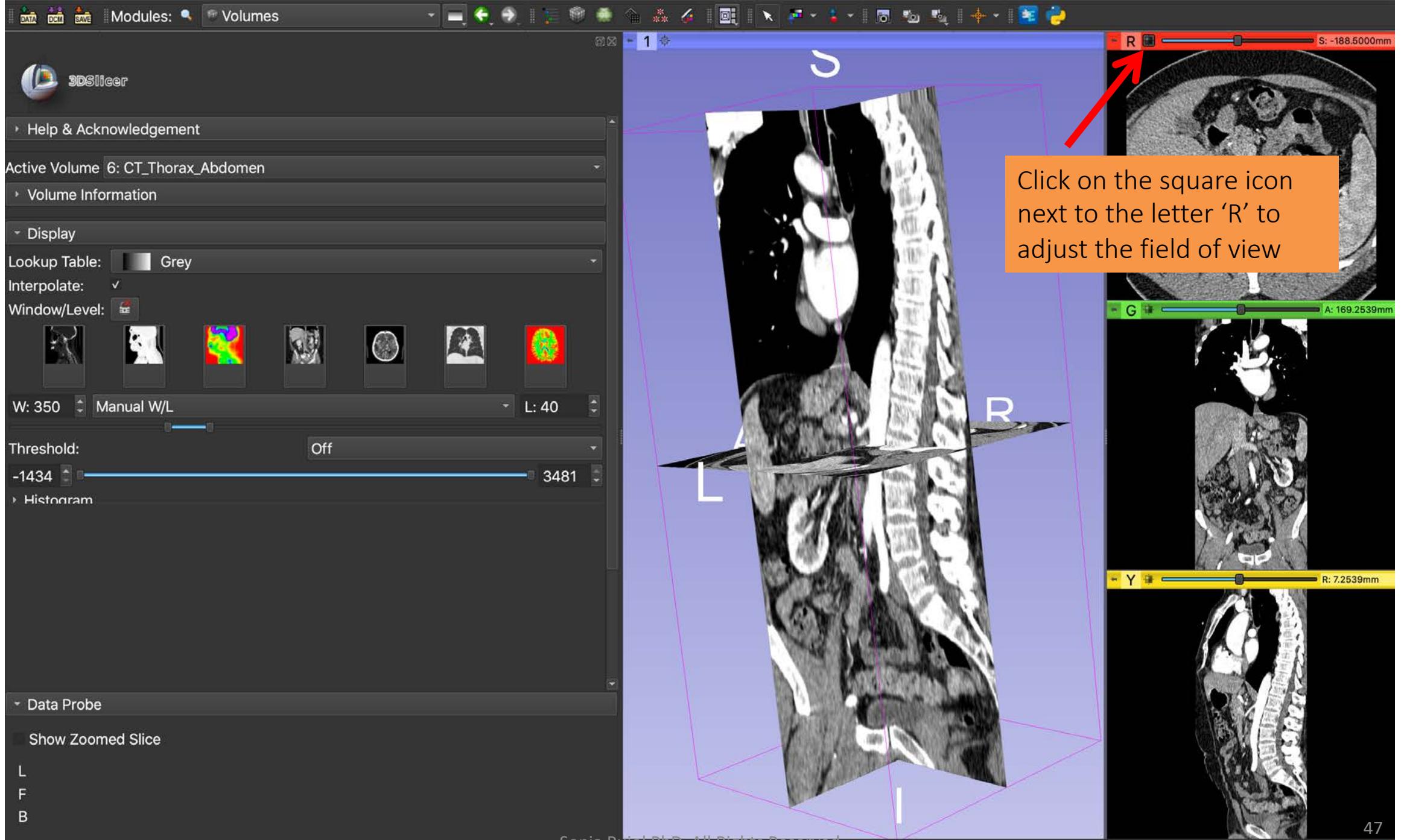

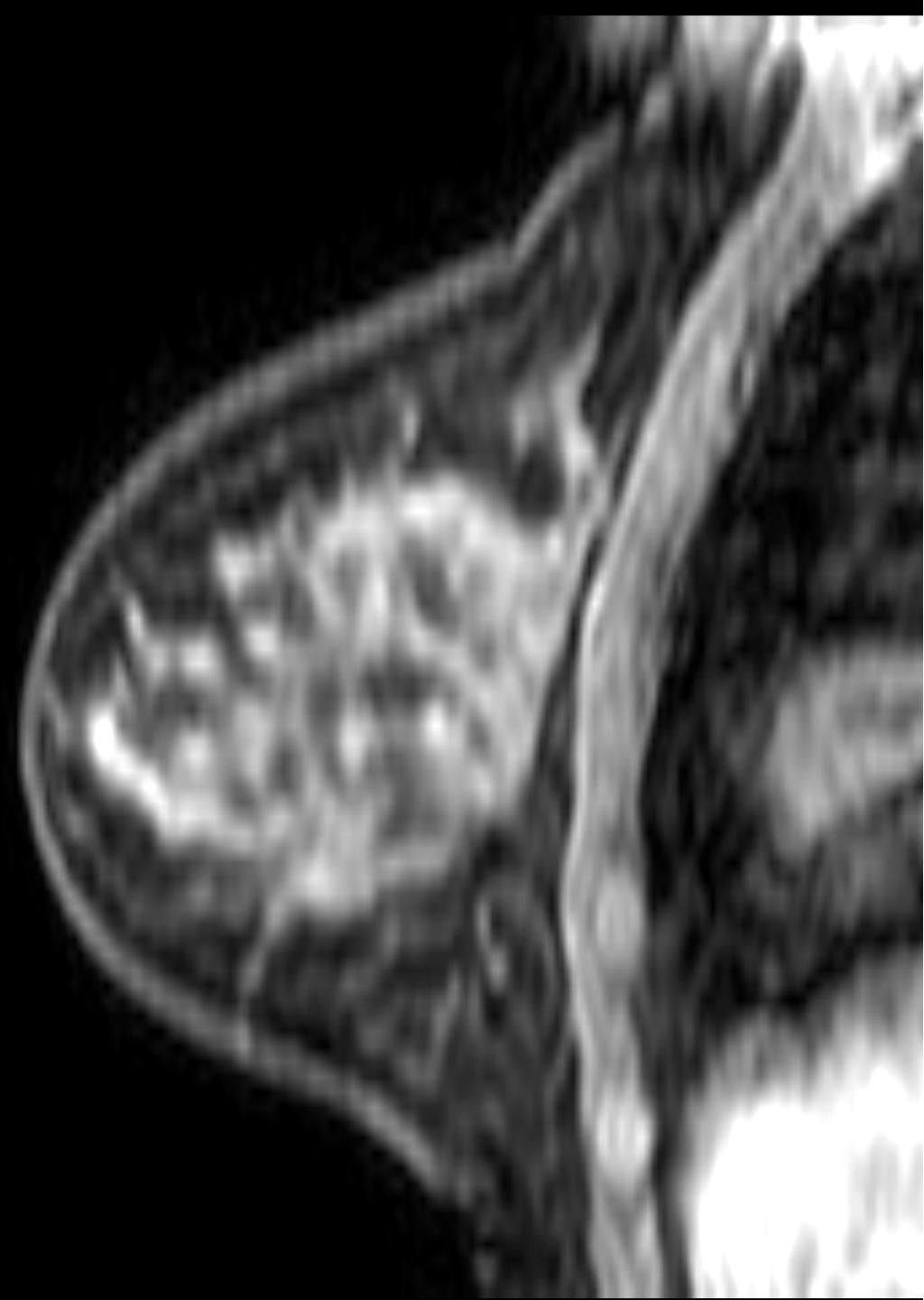


Dataset #2

Breast MRI

Breast MRI Dataset

- The Breast MRI dataset is part of the BREAST-DIAGNOSIS collection of The Cancer Imaging Archive (TCIA) of the National Cancer Institute
- The dataset was acquired on patient with right breast infiltrating ductal carcinoma
- The DICOM images consist of one study and three series: T2, STIR and BLISS

Bloch, B. Nicolas, Jain, Ashali, & Jaffe, C. Carl. (2015). Data From BREAST-DIAGNOSIS. The Cancer Imaging Archive. <http://doi.org/10.7937/K9/TCIA.2015.SDNRQXXR>

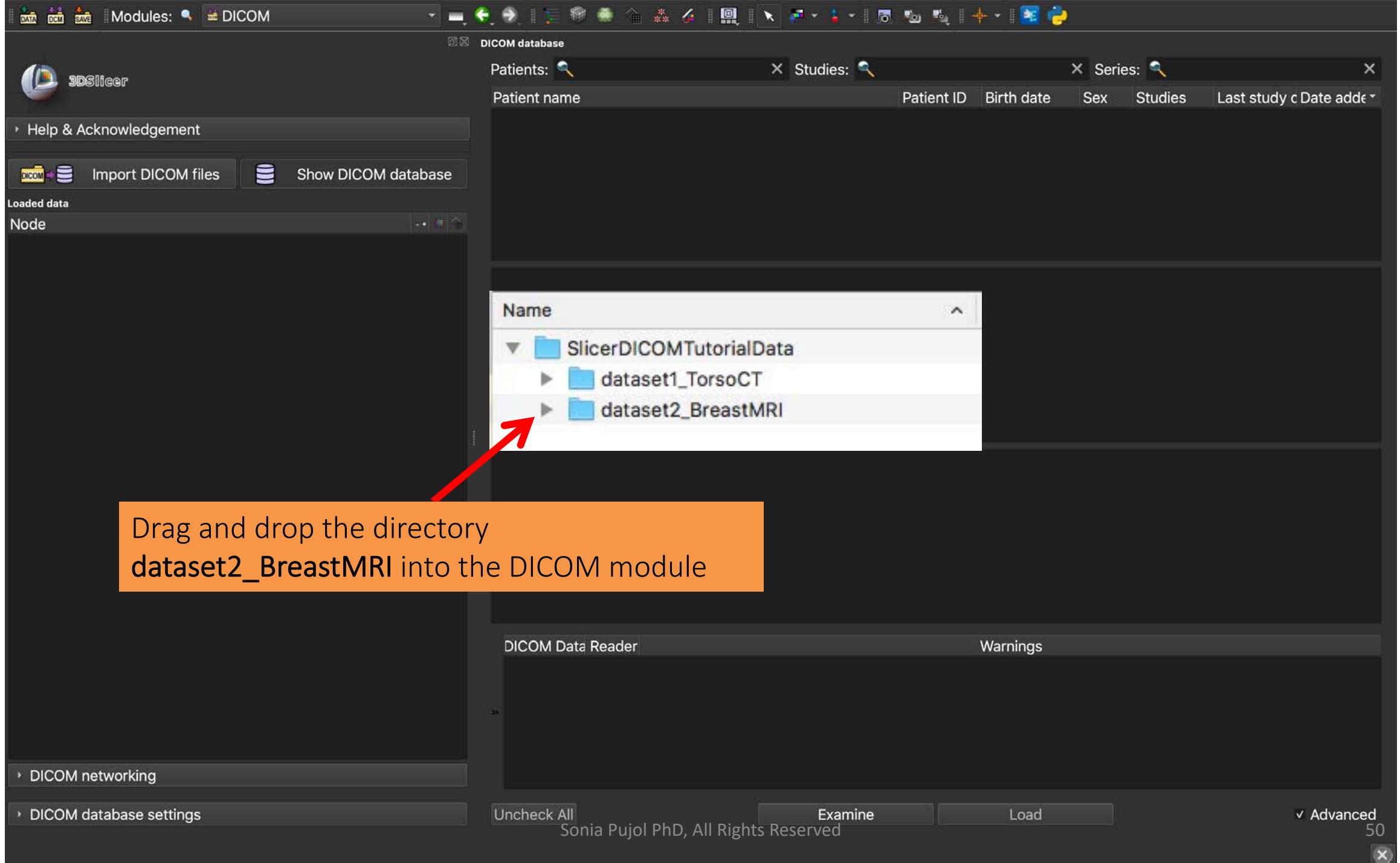

Slicer imports the **dataset2_BreastMRI** directory into the DICOM database

The directory contains 1 patient, 1 study and 3 series

 Import completed: added 1 patients, 1 studies, 3 series, 1008 instances.

OK

DICOM networking DICOM database settings

Uncheck All Examine Load Advanced

Sonia Pujol PhD, All Rights Reserved

51

Click on the PatientID
BreastDx-01-0005 to
display the study and the
three T2W, STIR and BLISS
series

Click on Examine

DICOM database

Patients: Studies: Series:

Patient name	Patient ID	Birth date	Sex	Studies	Last study	Date added
	BreastDx-01-0005		F	1	Tue Nov 11 2008	2020...583
patient1	patient1_ID			1	Wed Jun 1 2005	2020...273

Study date Study ID Study description Series Date added

20081111		MRI BREAST, BILATERAL WITH T WITHOUT CONTRAST	3	2020...583
----------	--	---	---	------------

Series # Series description Modality Size Count Date added

301	T2W_TSE SENSE	MR	528x528	84	2020...583
401	STIR SENSE	MR	528x528	84	2020...221
801	AX BLISS_AUTO SENSE	MR	528x528	840	202...646

DICOM Data Reader Warnings

DICOM networking

DICOM database settings

Uncheck All Examine Load Advanced

Sonia Pujol PhD, All Rights Reserved

52

DATA DCM SAVE Modules: DICOM

3DSlicer DICOM database

Help & Acknowledgement

Import DICOM files Show DICOM database

Loaded data

Node

- patient1 (patient1_ID)
 - CT Thorax Abdomen (20050601)
 - 6: CT_Thorax_Abdomen

Click on the double arrow to display the list of DICOM readers

DICOM database

Patient name	Patient ID	Birth date	Sex	Studies	Last study	Date added
	BreastDx-01-0005		F	1	Tue Nov 11 2008	2020...583
patient1	patient1_ID			1	Wed Jun 1 2005	2020...273

Study date	Study ID	Study description	Series	Date added
20081111		MRI BREAST, BILATERAL WITH T WITHOUT CONTRAST	3	2020...583

Series #	Series description	Modality	Size	Count	Date added
301	T2W_TSE SENSE	MR	528x528	84	2020...583
401	STIR SENSE	MR	528x528	84	2020...221
801	AX BLISS_AUTO SENSE	MR	528x528	840	202...646

DICOM Data	Reader	Warnings
301: T2W_TSE SENSE	Scalar Volume	
401: STIR SENSE	Scalar Volume	
801: AX BLISS_AUTO ...	MultiVolume	Images are not equally spaced (a difference of 2 vs 0 in spacings ...)
AX BLISS_AUTO SENSE	Scalar Volume	

Uncheck All Examine Load Advanced

Sonia Pujol PhD, All Rights Reserved

53

The screenshot shows the 3DSlicer application interface. At the top, there is a toolbar with various icons. Below the toolbar, the title bar reads "DICOM database". The main window is divided into several panels:

- DICOM database:** A table showing patient and study information. A red arrow points from the text "The list of DICOM plugins appear" to the list of plugins below this panel.
- Loaded data:** A tree view showing the hierarchy of loaded data, including "patient1 (patient1_ID)" and "CT Thorax Abdomen (20050601)".
- Node:** A panel showing the current node selected in the loaded data tree.
- Study description:** A table showing study details: Study date 20081111, Study ID MRI BREAST, BILATERAL WITH T WITHOUT CONTRAST, Series 3, Date added 2020-0..28.583.
- Series description:** A table showing series details: Series # 301, 401, 801, Series description T2W_TSE SENSE, STIR SENSE, AX BLISS_AUTO SENSE, Modality MR, Size 528x528, Count 84, Date added 2020-0..28.583, 2020-0..29.221, 2020-0..28.646.
- DICOM Data:** A table listing DICOM data entries with columns: DICOM Data, Reader, and Warnings. The entries are: 301: T2W_TSE SENSE (Scalar Volume), 401: STIR SENSE (Scalar Volume), AX BLISS_AUTO SENSE - as a 10 ... (MultiVolume), 801: AX BLISS_AUTO SENSE (Scalar Volume), AX BLISS_AUTO SENSE - as a 10 ... (MultiVolume), and AX BLISS_AUTO SENSE - as a 10 ... (MultiVolume). A warning message for the first entry states: "Images are not equally spaced (a difference of 2 vs 0 in spacing...)".

The list of DICOM plugins appear

- ✓ DICOMPETSVPlugin
- ✓ DICOParametricMapPlugin
- ✓ DICOMRWVMPPlugin
- ✓ DICOMScalarVolumePlugin
- ✓ DICOMSegmentationPlugin
- ✓ DICOMSlicerDataBundlePlugin
- ✓ DICOMTID1500Plugin
- ✓ DICOMVolumeSequencePlugin
- ✓ DicomRtImportExportPlugin
- ✓ DicomSrlmportExportPlugin
- ✓ MultiVolumeImporterPlugin

Uncheck All Examine Load Advanced

Sonia Pujol PhD, All Rights Reserved

Slicer DICOM Plugins

- ✓ DICOMPETSUVPlugin
- ✓ DICOMParametricMapPlugin
- ✓ DICOMRWVMPPlugin
- ✓ DICOMScalarVolumePlugin
- ✓ DICOMSegmentationPlugin
- ✓ DICOMSlicerDataBundlePlugin
- ✓ DICOMTID1500Plugin
- ✓ DICOMVolumeSequencePlugin
- ✓ DicomRtImportExportPlugin
- ✓ DicomSroImportExportPlugin
- ✓ MultiVolumeImporterPlugin

- Slicer implements a list of DICOM plugins to handle a diverse set of DICOM data objects
- These plugins need to be enabled in order to read specific DICOM data objects such DICOM RT data or DICOM DWI data

3DSlicer

Modules: DICOM

DICOM database

Patients: Studies: Series:

Patient name	Patient ID	Birth date	Sex	Studies	Last study date	Date added
patient1	BreastDx-01-0005		F	1	Tue Nov 11 2008	2020-0...28.583
	patient1_ID			1	Wed Jun 1 2005	2020-0...43.273

Import DICOM files Show DICOM database

Loaded data

Node

- patient1 (patient1_ID)
 - CT Thorax Abdomen (20050601)
 - 6: CT_Thorax_Abdomen

Study date Study ID Study description Series Date added

20081111		MRI BREAST, BILATERAL WITH T WITHOUT CONTRAST	3	2020-0...28.583
----------	--	---	---	-----------------

Series # Series description Modality Size Count Date added

301	T2W_TSE SENSE	MR	528x528	84	2020-0...28.583
401	STIR SENSE	MR	528x528	84	2020-0...29.221
801	AX BLISS_AUTO SENSE	MR	528x528	840	2020-0...28.646

✓ DICOMPETSUVPPlugin ✓ DICOParametricMapPlugin ✓ DICOMRWVMPPlugin ✓ DICOMScalarVolumePlugin ✓ DICOMSegmentationPlugin ✓ DICOMSlicerDataBundlePlugin ✓ DICOMTID1500Plugin ✓ DICOMVolumeSequencePlugin ✓ DicomRtImportExportPlugin ✓ DicomSrlImportExportPlugin ✓ MultiVolumeImporterPlugin


301: T2W_TSE SENSE ✓ 401: STIR SENSE ✓ AX BLISS_AUTO SENSE - as a 10 ... MultiVolume
801: AX BLISS_AUTO SENSE Scalar Volume Images are not equally spaced (a difference of 2 vs 0 in spacing...) AX BLISS_AUTO SENSE - as a 10 ... MultiVolume AX BLISS_AUTO SENSE - as a 10 ... MultiVolume



Uncheck All Examine Load Advanced

Sonia Pujol PhD, All Rights Reserved


Click on Load to load the data into Slicer


56

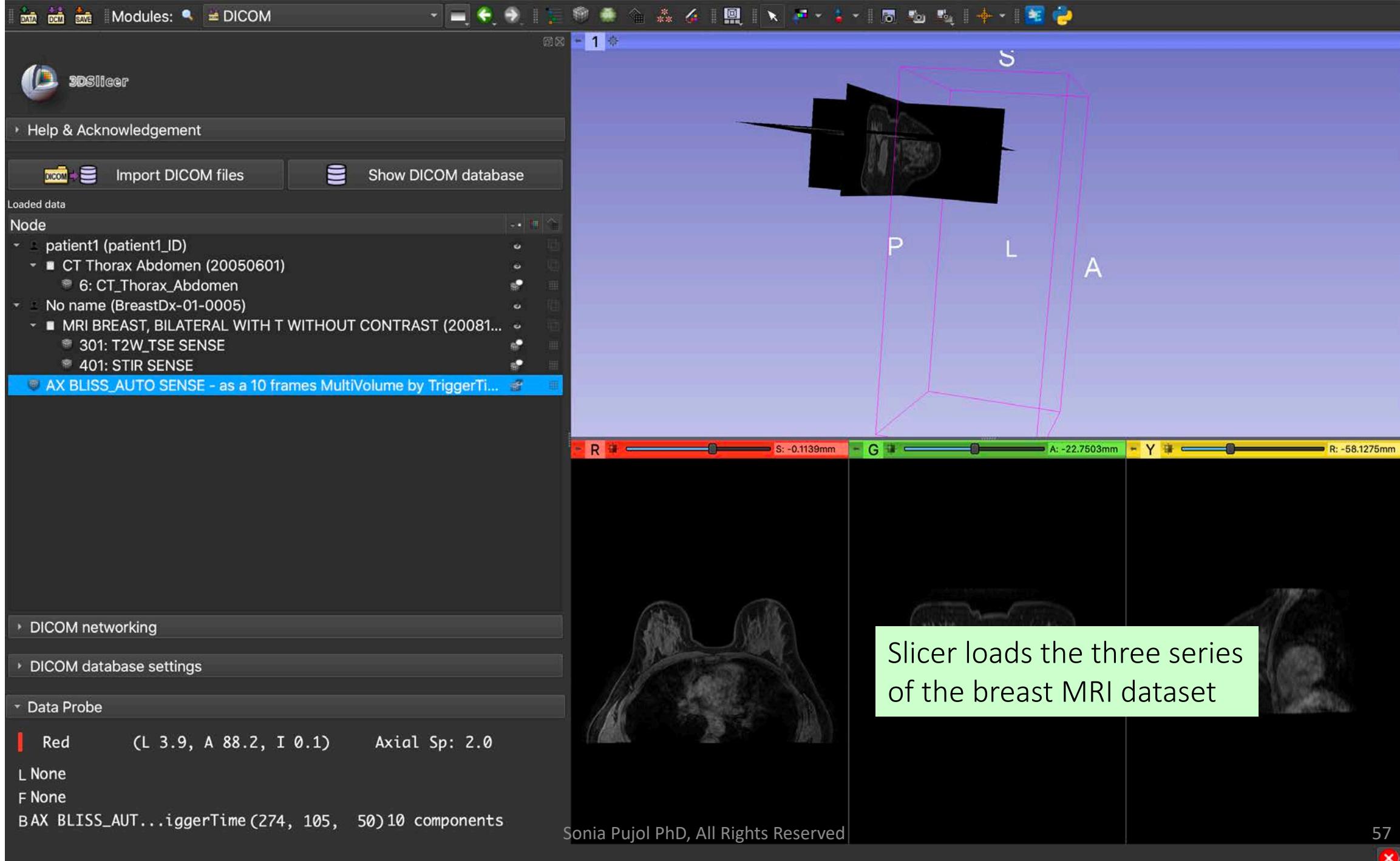

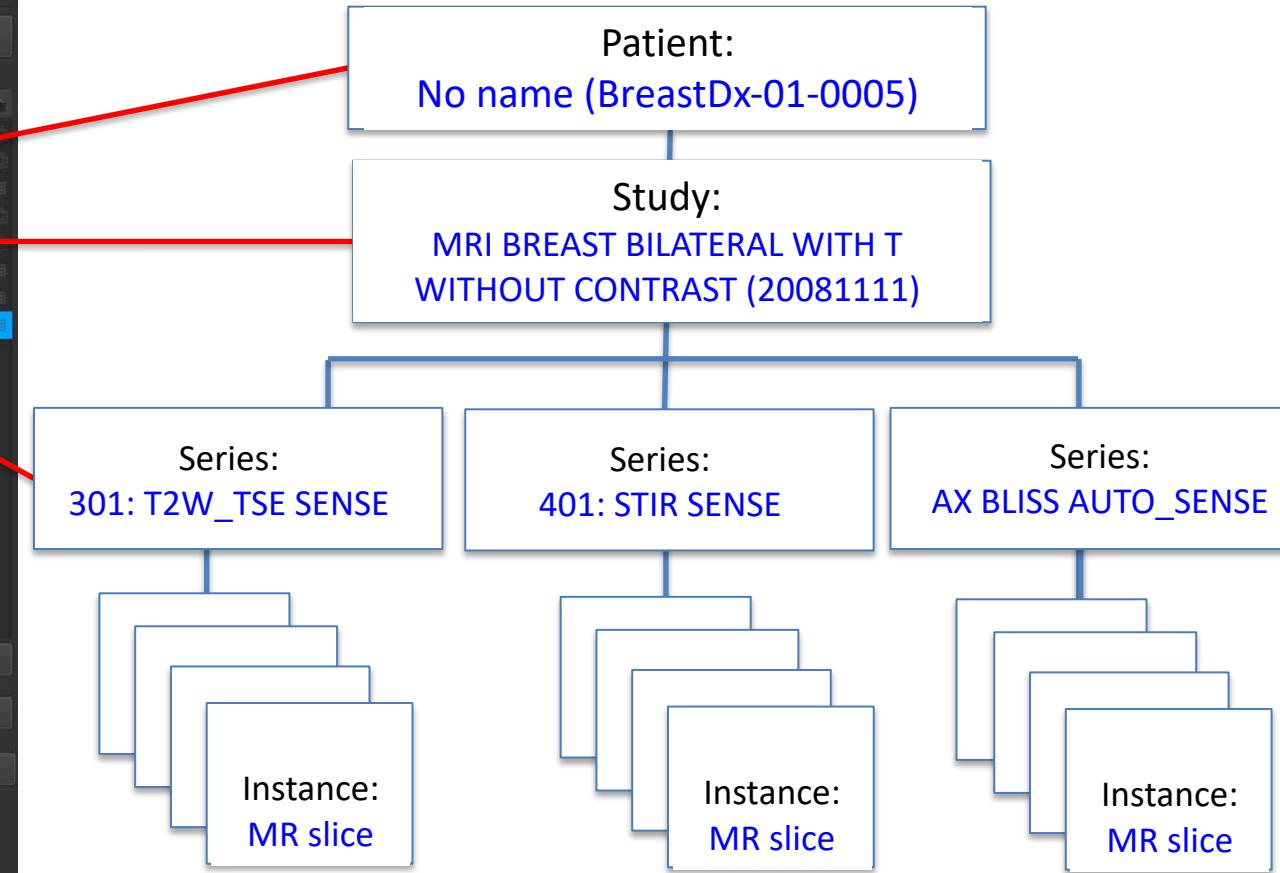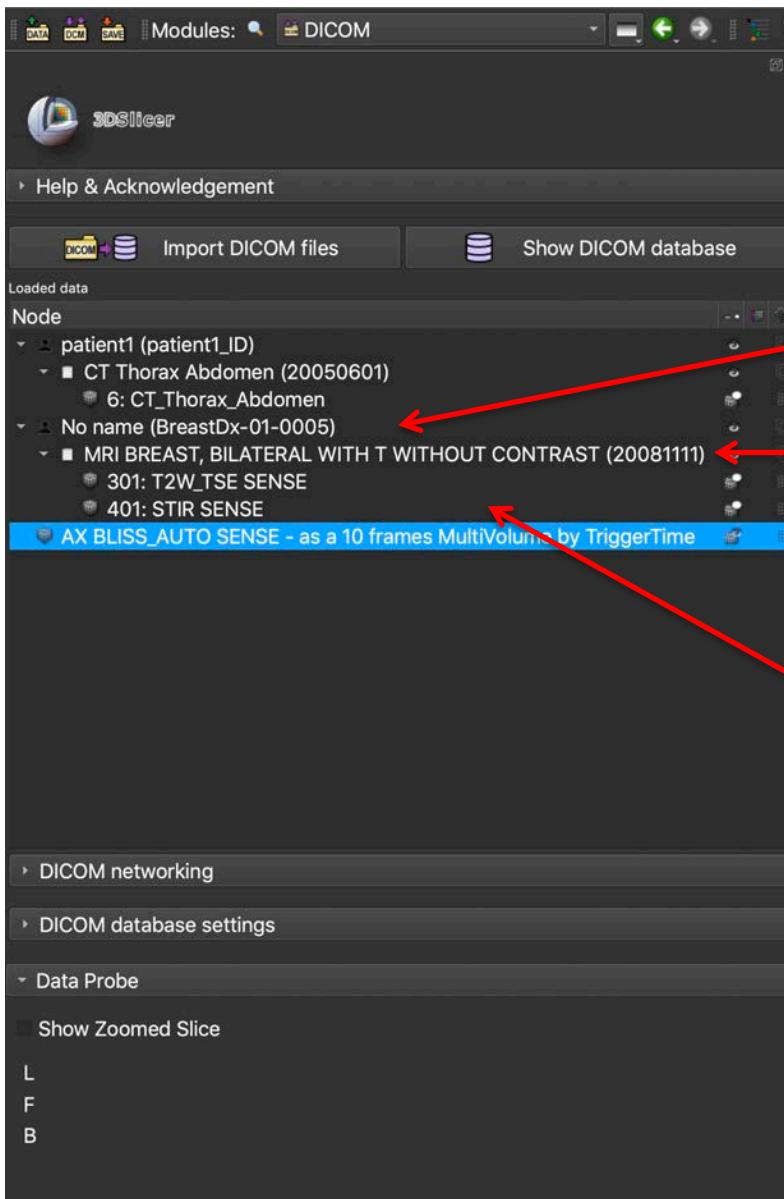

DICOM data are loaded into Slicer as a patient-study-series hierarchy

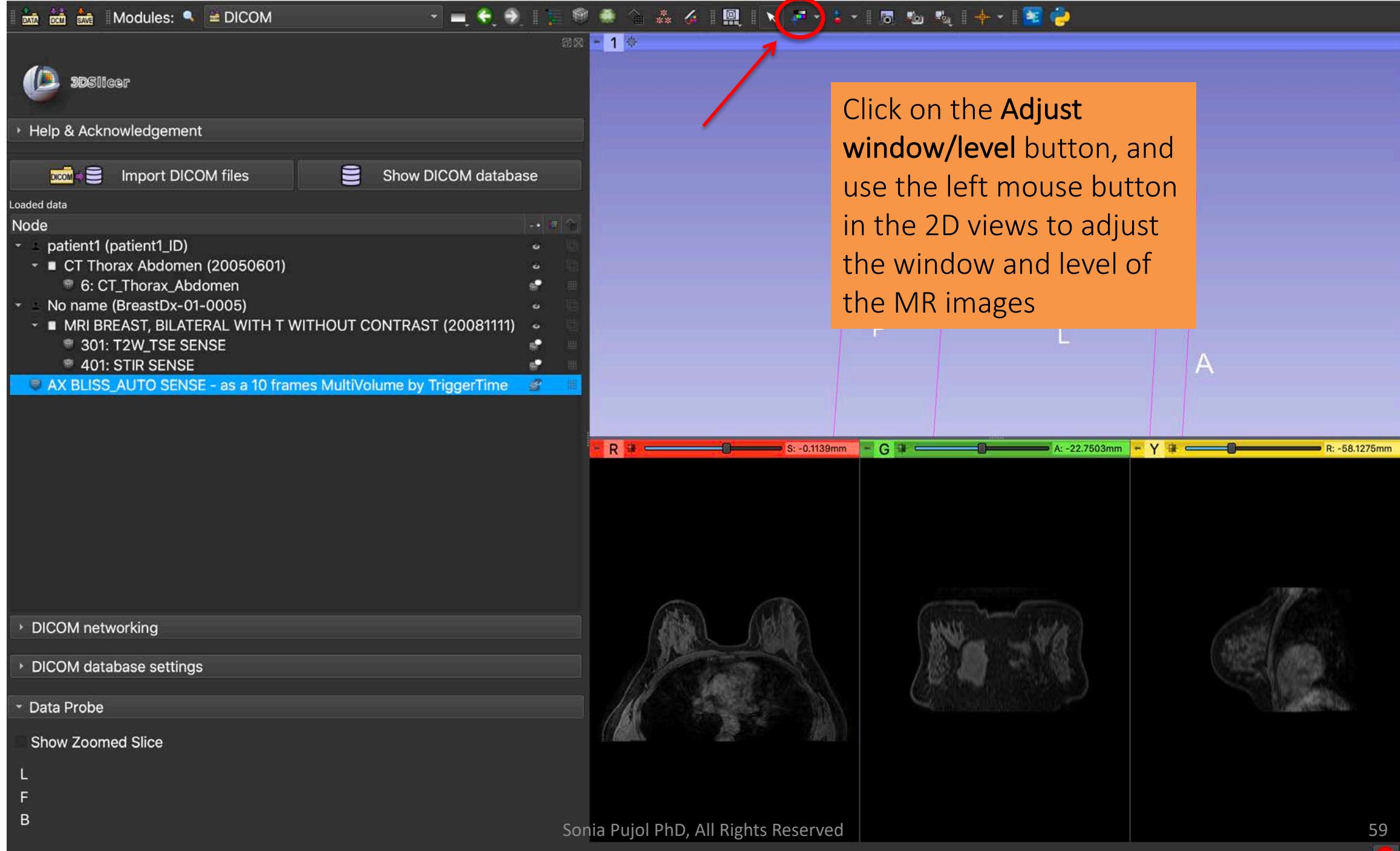

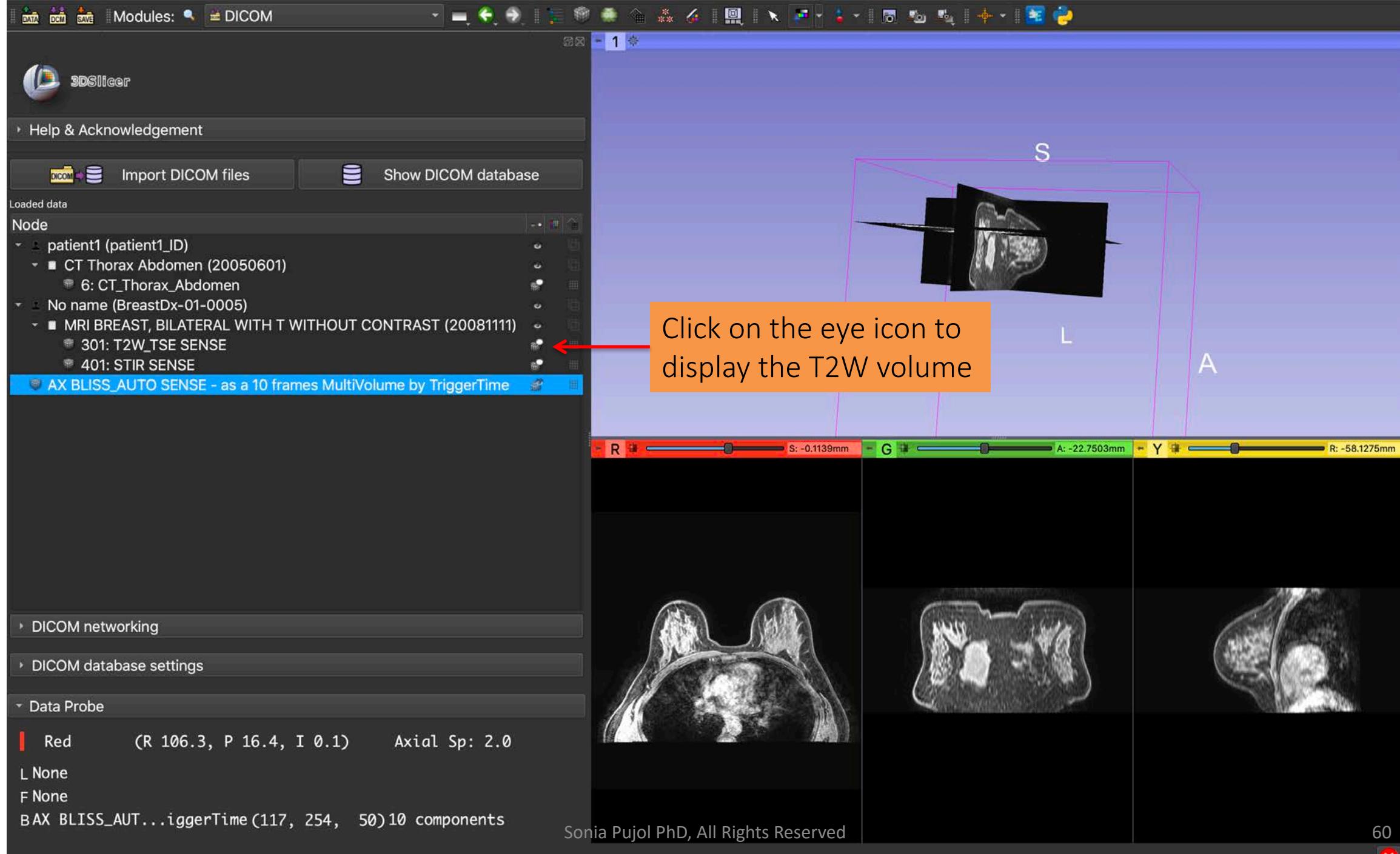



3DSlicer

Help & Acknowledgement


DICOM Import DICOM files

Show DICOM database

Loaded data

Node

- patient1 (patient1_ID)
 - CT Thorax Abdomen (20050601)
 - 6: CT_Thorax_Abdomen
 - No name (BreastDx-01-0005)
 - MRI BREAST, BILATERAL WITH T WITHOUT CONTRAST (20081111)
 - 301: T2W_TSE SENSE
 - 401: STIR SENSE
- AX BLISS_AUTO SENSE - as a 10 frames MultiVolume by TriggerTime

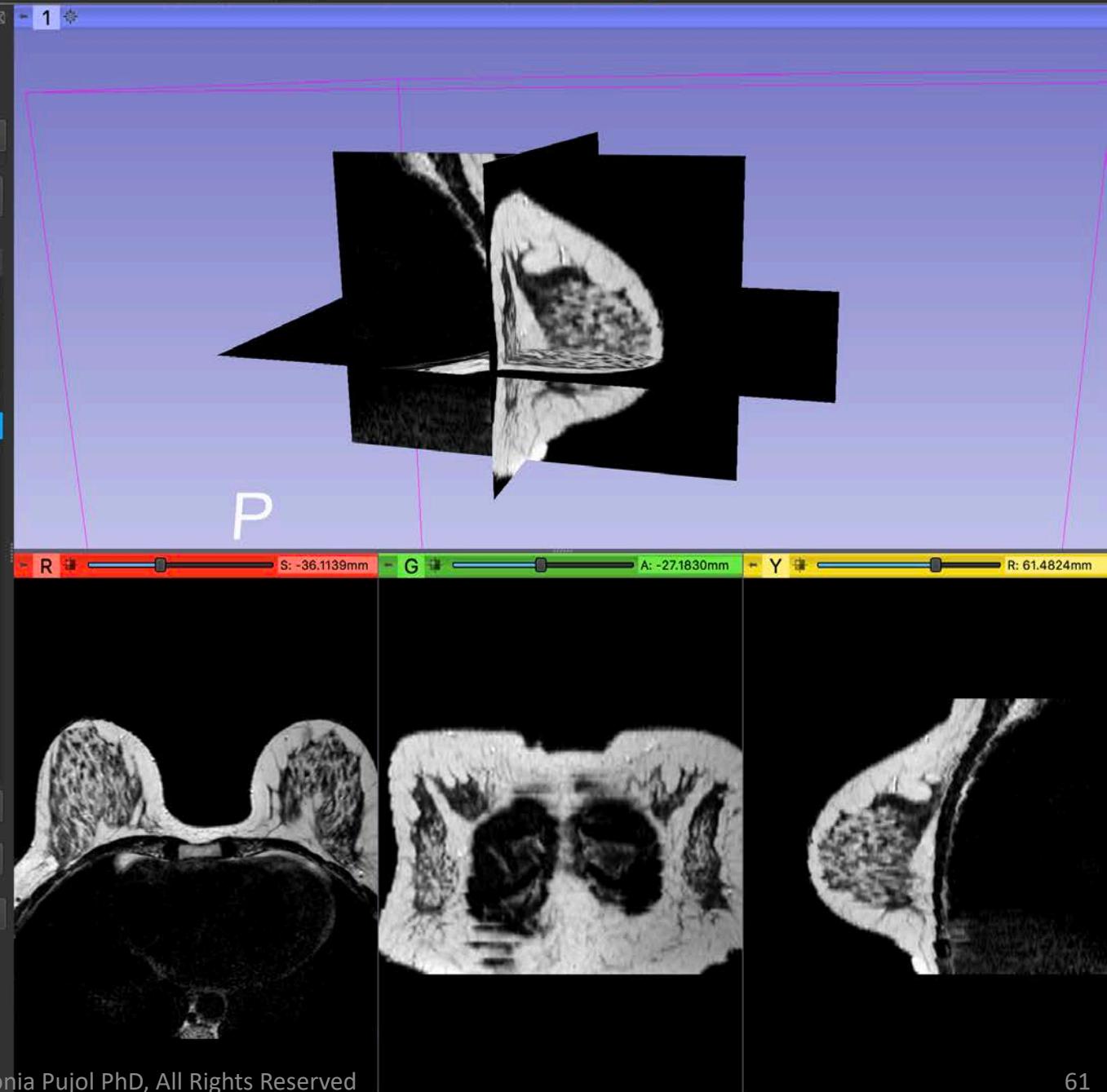

Slicer displays the
T2W_TSE volume in
the viewer

DICOM networking

DICOM database settings


Data Probe

Show Zoomed Slice

L

F

B

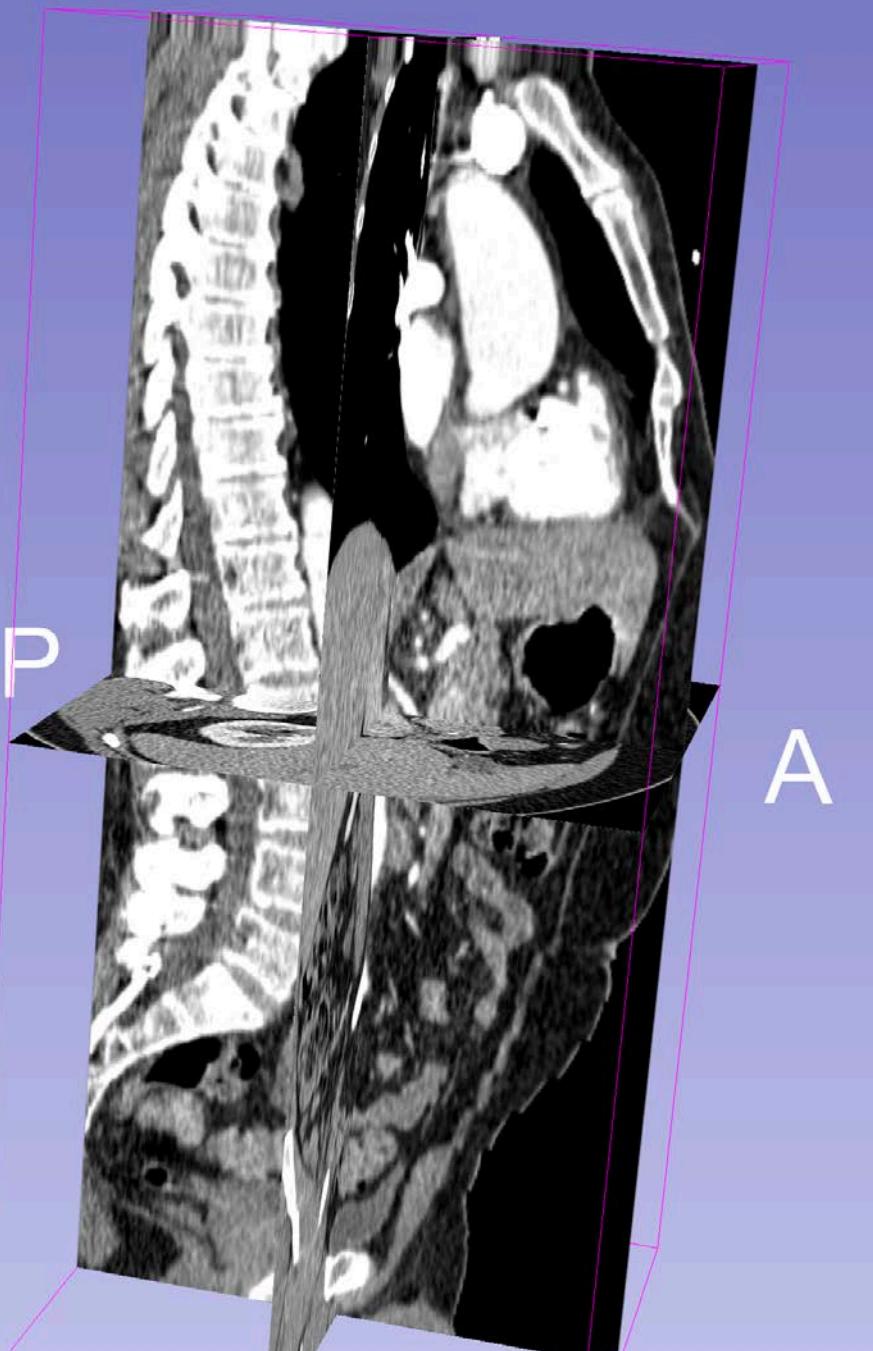

Conclusion

- 3D Slicer and the DICOM standard enable compliance with the FAIR principles for biomedical research
- By enabling interoperability between research and clinical environment, 3D Slicer and the DICOM standard lower the inherent barriers to the translation of research advances to patient care

Acknowledgments

 Neuroimage Analysis Center
(NIBIB P41 EB015902)